



# Cardiopulmonary Exercise Testing Special Edition

January 2019

# Editorial


## Editorial

In recent years cardiopulmonary exercise testing has become an increasingly more important tool and established itself as a valuable differential diagnosis in the fields of cardiology, respiratory and sports medicine. Furthermore, it is routinely used as a standard test method in other fields, such as industrial medicine, rehabilitation, and anaesthesiology for pre-operative risk assessment.

Due to advances in modern technology, administering a cardiopulmonary test has become easier, while the possibilities for evaluation and diagnosis have increased significantly.

Of special importance are individual measured values such as oxygen uptake as well as the graphic display of dynamic changes such as exercise flow-volume loops, the aerobic capacity, or the ventilatory efficiency ( $V'E/V'CO_2$ ) slope.

In this, our third special edition of cardiopulmonary exercise testing we would like to introduce our latest developments in this field as well as giving you further background information on topics such as: threshold determination and on testing with increased oxygen supplementation. We are pleased to contribute to the ongoing development of this important area of physiological measurement and we hope that this special edition will provide you with new and interesting insights into the world of cardiopulmonary exercise testing.

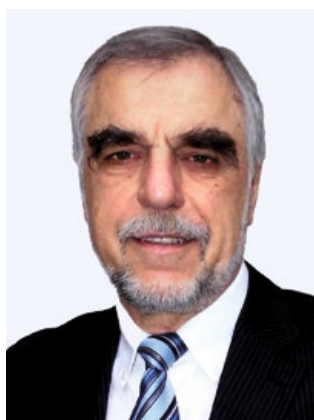

**Dr. Hermann Eschenbacher**  
Germany

## Table of Contents

### Editorial

Relevance of Cardiopulmonary Exercise Testing ..	3
Fields of Application .....	5
Evaluation .....	6

### Vyntus CPX - the Latest Product Generation

The Vyntus™ CPX at a Glance .....	12
Vyntus™ CPX - the Software .....	15

### Vyntus CPX - Options

Vyntus™ ECG .....	26
Vyntus™ CPX - High / Low FIO2 Option .....	28
Canopy - Indirect Calorimetry .....	30

### Basics and Diagnostics

Threshold Determination .....	32
Indirect Calorimetry .....	37
Haldane and Eschenbacher Transformation .....	40
Our CPET History .....	44

### Device Presentation

Vyntus™ WALK .....	46
--------------------	----

### The Last Page

Promotional Material .....	47
CPET Workshops .....	47

## Relevance of Cardiopulmonary Exercise Testing

### Introduction

**C**ardio **P**ulmonary **E**xercise **T**esting (CPET, sometimes abbreviated as CPX) is the determination of a person's performance during physical exercise by measuring, or calculating, the metabolic gas exchange alongside a number of other parameters.

In order to perform a specific task, the body needs to provide the required energy. This energy is primarily produced by the breakdown of carbohydrates, fats and proteins in the presence of oxygen (aerobic metabolism). So in addition to fuel (through food intake), the body needs to provide the muscles with sufficient oxygen for this metabolic process. The rate of oxygen required increases with the intensity of the exercise. As with any other burning process, carbon dioxide (CO<sub>2</sub>) is produced. This CO<sub>2</sub> is then transported from the muscle cells by the blood to the lungs, where it is removed from the body (via respiration).

When exceeding a certain level of exercise, the body will not be able to provide sufficient oxygen to all of the exercising muscles. The additional required energy is then produced by means of the so-called anaerobic metabolism. Due to the limits of the anaerobic energy reserves, the body will only be able to exercise at this level for a short period of time until it is exhausted or these reserves are depleted. Anaerobic metabolism also results in the production of additional CO<sub>2</sub> which is discarded through the lungs driving increased rates of ventilation.

Ultimately a number of finely tuned physiologic functions need to fit and interact together like a well-oiled machine, much like the impression the gear wheel model, created by Wasserman (simplified), illustrates so concisely below:

Ambient air is inhaled via the lungs and some of the oxygen present in the air (oxygen uptake  $\dot{V}O_2$ ) diffuses through the lung membrane into the blood where it is absorbed by hemoglobin and delivered to muscles by the cardiovascular system (via the circulation blood). Once in the muscle, the actual breakdown of substrate takes place, providing the patient with energy and enabling the body to perform mechanical work (exercise). The CO<sub>2</sub> produced during that process is also absorbed into the blood, transported back to the capillary blood of the lungs, across the lungs membrane, and finally exhaled (carbon dioxide production  $\dot{V}CO_2$ ).

Even from this simple description, it is easy to see that by measuring parameters such as ventilation,  $\dot{V}O_2$ ,  $\dot{V}CO_2$  and heart rate we can begin to determine the overall capacity of the system and start to pinpoint where any limitations may exist.

The aim of a standard CPET protocol is for the individual to be exposed to a load using a bike ergometer or a treadmill and incrementally increase workload for about 8 to 12 minutes until they can go no further. These are often referred to as an incremental ramp protocol to a volitional maximum. During the test, the patient is connected via mask (or mouth piece) giving minute ventilation, breathing frequency, oxygen uptake, carbon dioxide production as well as other parameters and heart rate is measured from an ECG of the measuring system. This procedure makes it possible to determine the maximum exercise capacity as well as various thresholds such as the endurance capacity threshold (also see chapter "Threshold Determination"). If an organ or organ-system is somehow impaired, the patient will fail to cope with the increasing load. In such a case, the characteristic patterns of the measured parameters can provide important information on which systems are affected by such impairment.



# Editorial

In sports medicine / science, step protocols are often used in order to receive more precise information regarding the speed or power of the athlete. When measuring patients, on the other hand, a ramp protocol is usually preferred as this will allow the patient to approach the maximum load within an acceptable time range before the need to terminate the test due to exhaustion (and not due to the maximum exercise capacity).

Thus, many aspects of cardiopulmonary exercise testing can be tailored to the individual's needs and capabilities.

With CPET it is possible to receive significant information on single functions or limitations. This leads to the key application areas for cardiopulmonary exercise testing:

- Determining the individual exercise capacity
- Determining the severity of a performance limitation
- Determining the aerobic and anaerobic performance ranges
- Determining and analyzing the effect of therapeutic interventions and / or rehabilitation in patients with performance limitation
- Differential diagnosis regarding possible causes for a performance limitation such as
  - Pulmonary limitation
  - Malfunction of the gas exchange
  - Cardiac limitation
  - Peripheral limitation
  - Motivational limitation

In contrast to simple stress tests, such as those performed with an ECG and a treadmill, CPET provides information on: test quality; allows for objective exercise capacity measurements; and points out causes for possible limitations. CPET also allows the possibility to assess the pre-operative risk for complications which may occur after a major surgery (such as lung or heart transplantation) more accurately. As a consequence, the post-surgical mortality rate can be reduced.

The various parameters produced by a CPET can be broadly categorized into the following types:

## Measurement and Stress parameters

- Tidal volume (VT)
- Breathing frequency (BF)
- Inspired or expired oxygen concentrations ( $FIO_2$ ,  $FEO_2$ )
- Inspired or expired  $CO_2$  concentrations ( $FICO_2$ ,  $FECO_2$ )
- Workload (Watt, respectively speed and elevation)
- Heart rate (HR, Stress ECG)
- Oxygen saturation ( $SpO_2$ )
- ...

## Calculated Parameters

- Respiratory minute ventilation ( $V'E$ )
- Oxygen uptake ( $V'O_2$ )
- Carbon dioxide production ( $V'CO_2$ )
- Respiratory Exchange Ratio (RER)
- Oxygen pulse ( $O_2$ -Pulse)
- Breathing equivalent ( $EqO_2$ ,  $EqCO_2$ )
- Dead space ventilation ( $VD / VT$ )
- Breathing reserve (BR)
- Heart Rate Reserve (HRR)
- ...

## Further evaluation parameters such as

- Threshold determination (VT1, VT2, VT3)  
(for further information, please see chapter "Threshold determination")
- Maximum oxygen uptake ( $V'O_{2max}$ )
- Slope determination
- Aerobic capacity ( $dV'O_2 / dWR$ )
- Ventilatory efficiency ( $V'E(V'CO_2)$  slope)
- Alveolar-arterial oxygen pressure difference ( $P(A-a)O_2$ )
- ...

## Fields of Application

The various parameters, measured and calculated, turn CPET into a comprehensive and highly informative method with applications in numerous fields of medicine:

### *Respiratory Medicine*

- Obstructive and restrictive ventilatory disorders
- Interstitial disorders
- Pulmonary hypertension
- Diffusion and distribution disorders
- Flow limitations
- Exercise related dyspnea of unknown origin
- Suspected limited exercise capacity due to circulatory or pulmonary vascular disorders
- Suspected exercise-induced asthma
- Trending for subtle respiratory disease changes
- Pre-operative risk assessment for lung transplant patients

### *Cardiology*

- Coronary heart disease
- Cardiomyopathy
- Heart disease, valvular heart failure
- Congenital cardiac defects
- Pre-operative risk assessment for heart transplant patients
- Cardiac insufficiency

### *Sports Medicine / Science*

- Measurement of physical exercise capacity
- Threshold determination
- Training management
- Quantification of training success

### *Occupational Medicine*

- Exercise-related career proficiency tests
- Determining the degree of disability or work limitation / inability
- Fitness checkups (high altitude, air travel, tropical climate, diving)

### *Intensive Care*


- Pre-operative risk assessment
- Nutrition control (adjusting parenteral nutrition of intensive care patients)

### *Rehabilitation*

- Optimizing rehabilitative measures
- Assessing and documenting rehabilitative and therapeutic progress

### *Nutrition*

- Determination of Resting Energy Expenditure
- Energy Expenditure during Exercise
- Substrate utilization
- Nutritional counselling
- Dietary advice

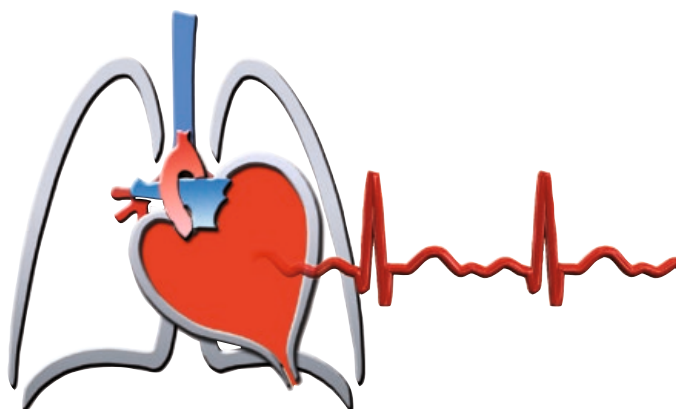


## CPET Evaluation

Different evaluation and interpretation procedures are used depending on whether the subject is a healthy athlete or a patient with cardiac and/or pulmonary limitation. The following considerations cannot be considered comprehensive but are intended to describe only the main aspects of a CPET evaluation. For more detailed information, please refer to the list of additional literature at the end of this edition.

When using CPET equipment it is desirable to be able to examine different parameters and graphs at different times, both during and after the measurement. The presentation of the data should be concise, comprehensive and systematic. The layout below uses the internationally recognized 9-Panel-plot according to Prof. Karlman Wasserman (Wasserman 2009). In 2012, the 9-Panel-plot was updated (Wasserman (2012)). Some of the panels have been moved around, but the information content remains the same. The following 9-panel-graphic considerations refer to the original order. The software, of course, allows the user to choose between both alternatives, or create their very own 9-panel layout.

The 9-Panel-plot offers a concise overview of the cardiovascular, ventilatory and gas exchange parameters:

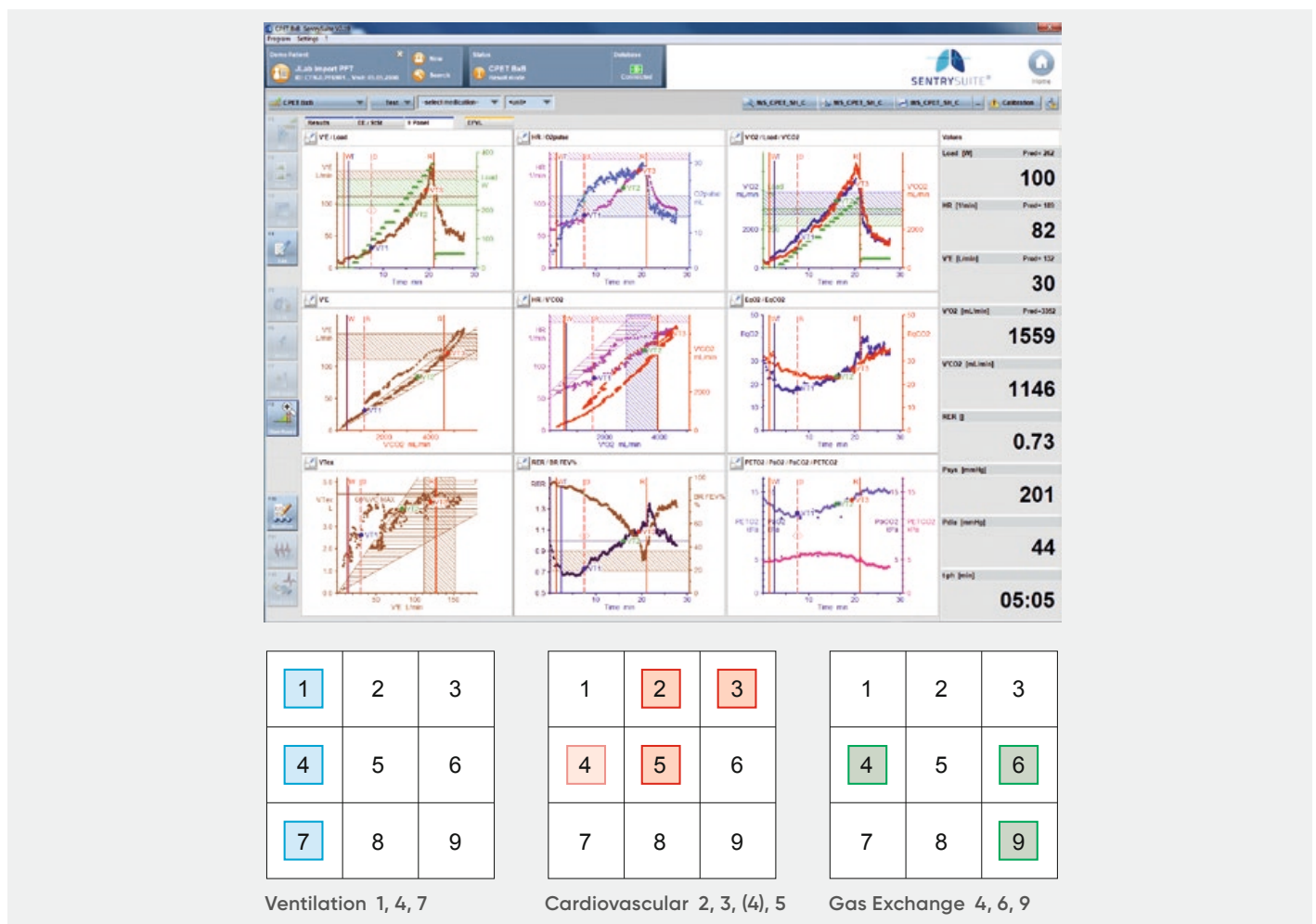

The blue panels mainly illustrate ventilatory aspects, the red ones relate to the cardiovascular parameters and the green panels convey gas exchange information.

Panel 1: V'E and load against time

Panel 2: HR and O<sub>2</sub> pulse against time

Panel 3: V'O<sub>2</sub>, V'CO<sub>2</sub> and load against time

Panel 4: V'E against V'CO<sub>2</sub>

Panel 5: HR and V'CO<sub>2</sub> against V'O<sub>2</sub>

Panel 6: EqO<sub>2</sub> and EqCO<sub>2</sub> against time

Panel 7: VTex against V'E


Panel 8: RER and BR FEV% against time

Panel 9: PETO<sub>2</sub> and PETCO<sub>2</sub> as well as PaO<sub>2</sub> and PaCO<sub>2</sub> against time

Cardiology Aspects

Panel 3


This panel illustrates the patient's general exercise performance and gives insight into oxygen delivery to, and utilization at, the exercising muscles. From this panel it is immediately apparent whether the subject has reached, or even exceeded, their expected exercise capacity (indicated by the hatched areas). If the expected exercise capacity is reached, it is safe to exclude a severe limitation. The oxygen uptake (in blue) of a healthy individual increases linearly with the workload (in green) in an approximate ratio of 10 mL/W. Providing the vertical axes are scaled to the same ratio (200 Watt = 2000 mL/min  $\dot{V}O_2$  or 200 Watt = 2.0 L/min  $\dot{V}O_2$ ) workload and  $\dot{V}O_2$  should increase in parallel to one another, though  $\dot{V}O_2$  will sometimes flatten upon reaching the peak work load. If flattening of the oxygen uptake occurs before the estimated exercise capacity (hatched area) is reached it is likely as a result of poor oxygen delivery to the muscles and evidence of cardiovascular limitation. It is not unusual for the  $\dot{V}O_2$  to flatten in athletes but this will occur at a level well above the expected  $\dot{V}O_2$  peak for a normal person. Furthermore, the panel provides information as to whether the peripheral muscle cells are utilizing sufficient oxygen. If this is not the case, the oxygen uptake will not increase linearly with increasing work load and will show a lower slope (less than 10 mL/W).

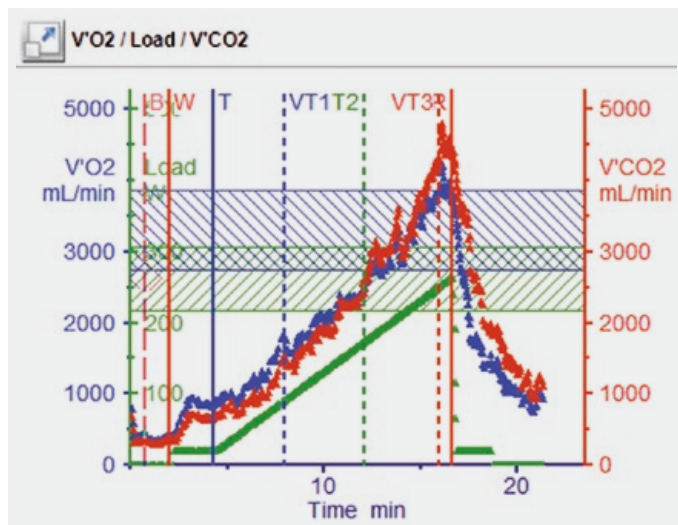

Panel 3 of the 9-Panel-Graphics. The respective phases are marked by vertical lines, the dashed lines indicate the respective thresholds and the hatched areas indicate the predicted values to be reached.

Panel 2


This panel reveals information on the patient's heart rate (HR) and oxygen pulse ( $O_2$  pulse). In healthy subjects, the heart rate is expected to rise with the increasing work load and will show a slight decline of the slope after some time, whereas patients with a cardiac impairment usually show a larger increase of the heart rate. With good cardiac function, the amount of oxygen transported per heart beat ( $O_2$  pulse) is high and increases throughout the test. Patients with poor cardiac function, the oxygen transport can only be increased by additional oxygen extraction. The oxygen pulse will reach a plateau as soon as this maximal extraction is reached. Consequently, a further increase in work load will result in a disproportionate increase of the heart rate.

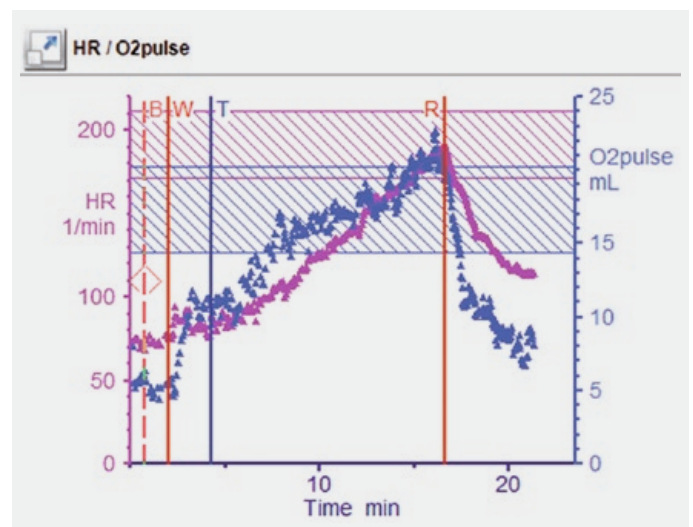

Panel 2 of the 9-Panel-Graphics.

Panel 5

In healthy subjects, the heart rate against oxygen uptake trace (in pink) will increase linearly as illustrated in this panel. In general the main areas of aerobic conditioning (increased stroke volume; higher mitochondrial density in the exercising muscles; and increased capillarization of those muscle) result in increased oxygen delivery and utilization and as a result the heart has to beat less to deliver oxygen, this is seen as a lower HR versus  $\dot{V}O_2$  slope. In deconditioned subjects (e.g. with low stroke volume) the converse is true and if there is an acute cardiac impairment this will be reflected by a sudden, disproportionate increase in heart rate. Panel 5 also displays the  $\dot{V}CO_2$  against  $\dot{V}O_2$  slope allowing the user to determine the different thresholds via the 'V-Slope' method (Beaver (1986)).

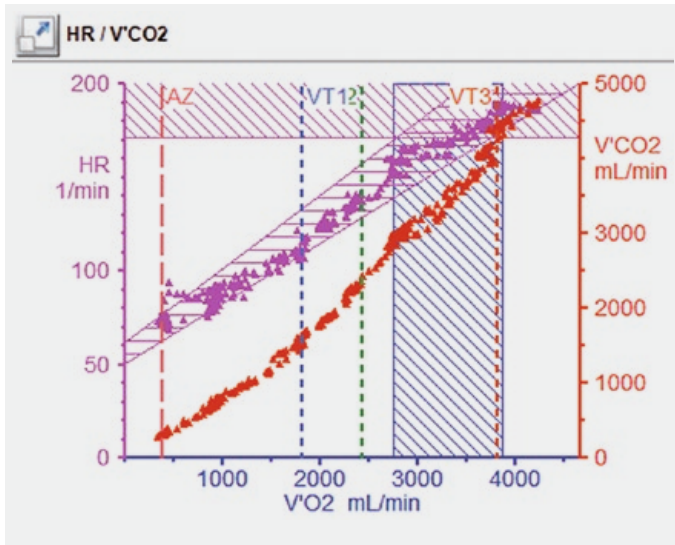

Panel 5 of the 9-Panel-Graphics

## Panel 4

Panel 4 provides cardiovascular information regarding, in particular, the pulmonary vascular circulation and for that reason will be described in detail both in the pulmonary as well as in the gas exchange section.

## Pulmonary Aspects


With panels 1, 4 and 7 it is possible to assess the ventilatory performance. In order to determine the maximal ventilation, both tabular and individual predicted values are of significance.

## Panel 1

This plot presents minute ventilation ( $V'E$ ) and workload (Watts) against time. In healthy subjects, the ventilation initially increases in a linear fashion. As exercise continues the trace increases out of proportion as it passes the respective ventilatory thresholds, this is caused by the increase in anaerobically produced  $CO_2$  and the ensuing metabolic acidosis (see "Threshold Determination"). This will only occur provided there is sufficient breathing reserve to accommodate this hyperpnea. In subjects suffering from pulmonary disease it is useful to display the subject's maximum ventilation obtained by means of a forced spirometry measurement (usually calculated from  $35 \times FEV1$ ) or a maximal voluntary ventilation maneuver (MVV) in order to detect ventilation limitation.


## Panel 4

This panel demonstrates the relationship between minute ventilation ( $V'E$ ) and the carbon dioxide production ( $V'CO_2$ ). A healthy subject requires an increase in ventilation ( $V'E$ ) of

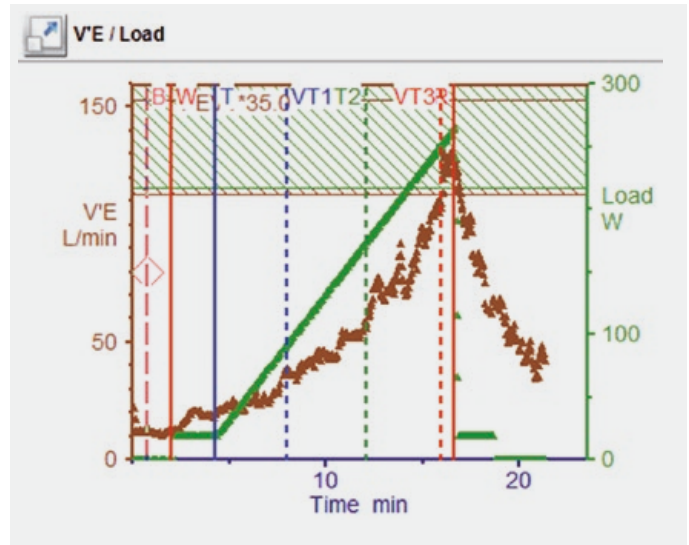

Panel 1 of the 9-Panel-Graphics. The respective phases are marked by the vertical lines, the dashed lines indicate the respective thresholds and the hatched areas indicate the predicted values to be reached.

about 25 L per additional liter of  $CO_2$ . If dead space ventilation is increased and/or an impairment of the gas exchange is present, the ventilation must be increased in order to expel the same amount of  $CO_2$ . Increased dead space ventilation shifts this curve upwards without increasing the slope, while an impaired diffusion results in a steeper slope. In this panel, it is also possible to display the maximum ventilation ( $35 \times FEV1$ ) next to the predicted values in order to easier determine whether there is any breathing reserve (BR).  $V'E$  and  $V'CO_2$  are closely tied and as a result this relationship is highly linear for much of the test, but when it reaches the ventilatory thresholds  $VT2$  and  $VT3$  (provided  $VT3$  has been reached) the slope increases due to the consequent hyperventilation.

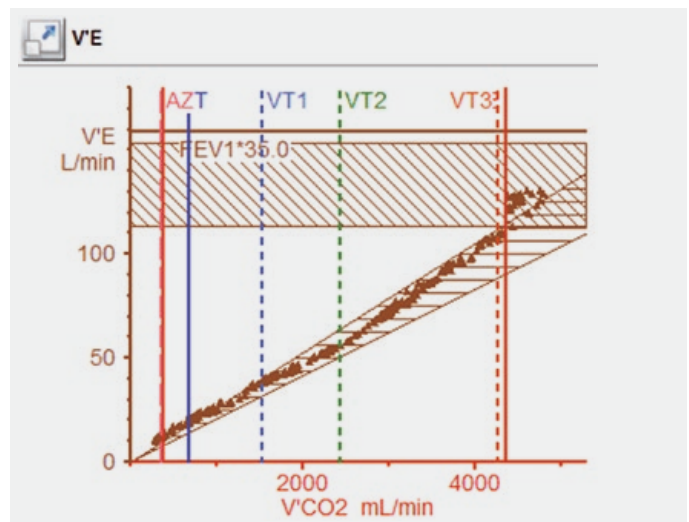


Panel 4 of the 9-Panel-Graphics.



**Panel 7**

This panel traces the changes in breathing pattern by plotting the expiratory tidal volume (V<sub>Tex</sub>) against the minute ventilation (V<sub>E</sub>). Unusual values suggest the presence of an obstructive or restrictive disorder. Patients with flow limitations will try to breathe as deeply and as slowly as possible which will cause the trace to curve along the upper isopleth (the straight line from the origin, in this case representing a breathing frequency of 20 breaths per minute). If a restrictive lung disorder is present, the patient will quickly reach the maximum tidal volume due to the low vital capacity. Further increases in ventilation are through possible only by increasing the breathing frequency. As a result, the curve will reach an early plateau and then run horizontally to intersect the lower isopleth (representing 50 breaths per minute). In addition to the predicted value of V<sub>E</sub>, the MVV value (maximum voluntary ventilation) and/or the patient's predicted value calculated from FEV<sub>1</sub>\*35 can be displayed in this panel to illustrate whether the patient has reached the maximum ventilation and whether a ventilatory impairment is present or not. Displaying the inspiratory capacity (IC) can also be useful as it corresponds to the maximum attainable tidal volume during the exercise. If the IC value has not been determined it can be approximated as 60% of the subject's vital capacity.

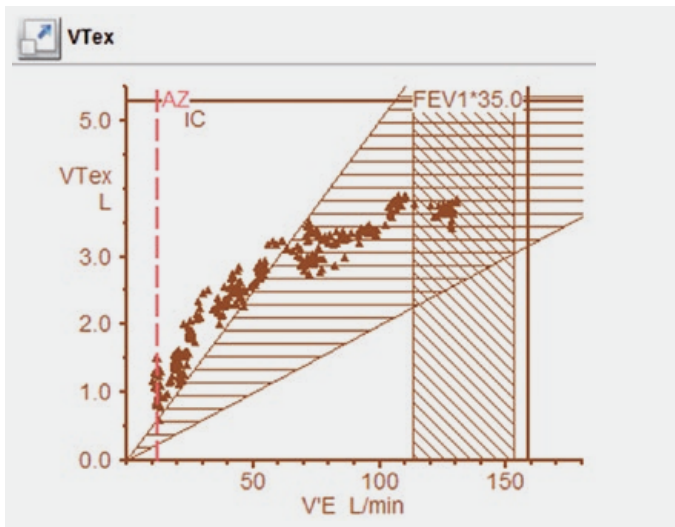


Panel 7 of the 9-Panel-Graphics.

**EFVL Measurement**


Another pulmonary aspect is dynamic hyperinflation which can be clearly demonstrated by measuring the flow-volume-curve during exercise (Exercise Flow Volume Loop - EFVL). This feature was already implemented into our previous version of software at the beginning of 1990. However, it is unfortunately not (yet?) considered in the 9-Panel layout.

By means of the EFVL measurement it is possible to immediately recognize a potential flow limitation of the subject which is indicated by the curve measured during exercise (blue) approaching or even slightly exceeding the maximum F/V curve (black) obtained from resting spirometry.

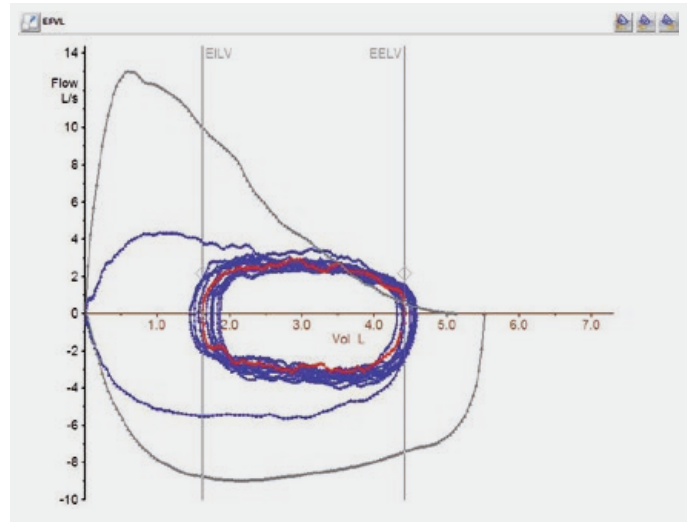

F/V curve during exercise (blue) compared to the maximum F/V curve (black)

If the EFVL curve is recorded several times during exercise, it is very easy to read if dynamic hyperinflation is present because the End-Expiratory Lung Volume (EELV) increases and/or if the subject exhausted his/her maximum tidal volume (End Inspiratory Lung Volume - EILV increases almost to TLC).

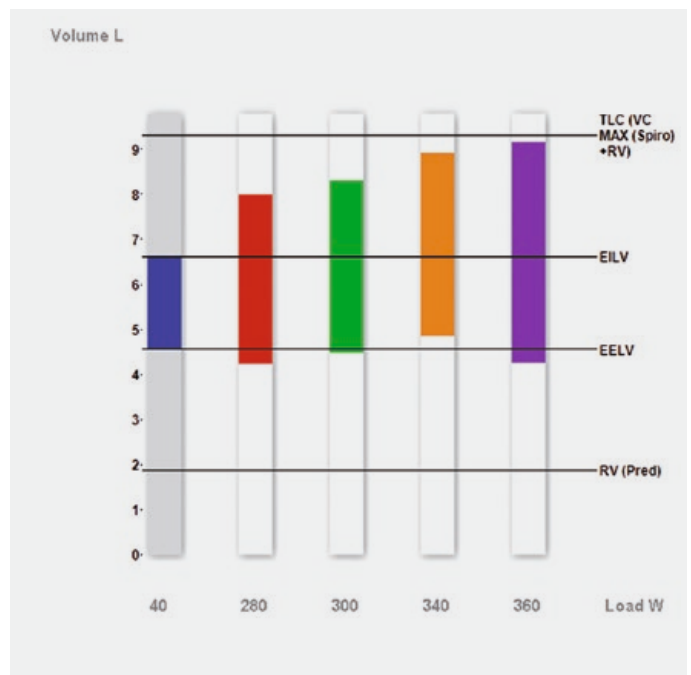


Exercise Flow Volume Loops recorded during CPET displayed as a bar diagram

# Editorial

## Gas Exchange

### Panel 6

Both panel 4 as well as panel 6 provide important information on the gas exchange. Panel 6 displays the breathing equivalents for  $\dot{V}O_2$  and  $\dot{V}CO_2$  ( $EqO_2$  and  $EqCO_2$ ). Please note, although the equivalents are approximately the same,  $EqO_2$  does not equal  $\dot{V}E/\dot{V}O_2$  and  $EqCO_2$  does not equal  $\dot{V}E/\dot{V}CO_2$ . This is because the breathing equivalents need to be corrected for the apparatus dead space, but despite this they are often incorrectly represented as  $\dot{V}E/\dot{V}O_2$  or  $\dot{V}E/\dot{V}CO_2$  in many publications. They give a measure of instantaneous ventilatory and gas exchange efficiency: How many L does the respective patient have to breathe in order to uptake 1 L oxygen or to produce 1 L carbon dioxide? At the beginning of the measurement, the values are relatively high due to the high dead space to tidal volume ratio ( $V_T$  low) and will decrease with the load as the tidal volume increases.  $EqO_2$  will reach a minimum (indicating optimum efficiency) in the  $VT1$  area,  $EqCO_2$  between  $VT1$  and  $VT2$ . Because of this, panel 6 can be used to help determining the ventilatory thresholds. A healthy person has a ventilatory demand of approximately 20–25 L in order to absorb 1 L oxygen and needs to ventilate approximately 25–30 L to release 1 L carbon dioxide. Elevated values indicate an inefficient gas exchange which can be caused by both an increased dead space ventilation and/or an impaired gas diffusion.

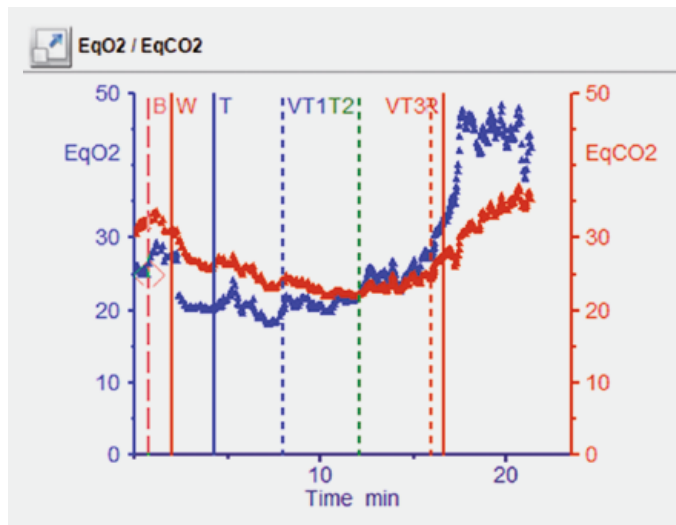


**Panel 6 of the 9-Panel-Graphics.** The respective phases are marked by vertical lines, the dashed lines indicate the respective thresholds.

### Panel 4

The respiratory minute ventilation  $\dot{V}E$  usually increases linearly until reaching  $VT2$  as the respiratory drive is primarily determined by the production of  $CO_2$  (the more  $CO_2$  is

released the higher the respiratory minute volume has to be). It corresponds (though not exactly) to the minimum  $EqCO_2$  and is approximately 25 L/L for  $CO_2$  in a healthy subject. Slopes of more than 40 L/L indicate additional cardiac impairment like pulmonary hypertension, thus this slope offers further valuable data for cardiologists.

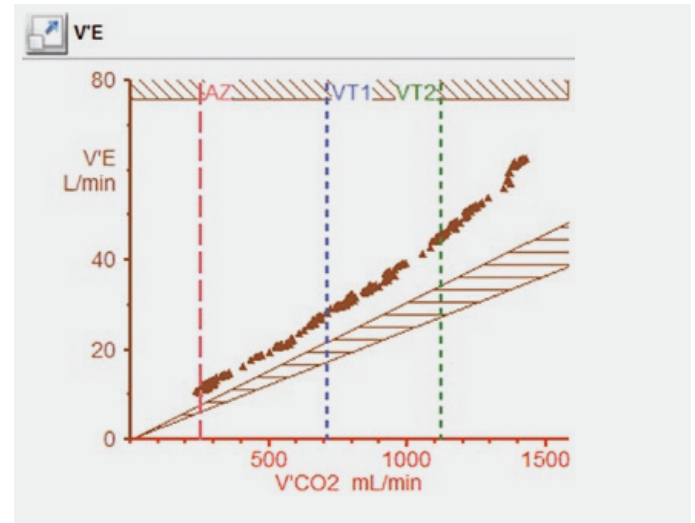

**Panel 4 of the 9-Panel-Graphics.** The respective thresholds are marked by the vertical lines. The inclining hatched area indicates the normal slope


In addition to the slope, a shift of this curve is important as well: With an increased dead space ventilation the patient needs to breathe more frequently from the start. Since the dead space usually does not change with increasing ventilation (but rather decreases due to the higher breathing volume) an increased dead space ventilation is indicated by an upwards shift of the curve.

Consequently, increased slopes indicate a diffusion impairment whereas a shift upwards is due to an increased dead space ventilation.

### Panel 9

Further information on gas exchange as well as on threshold determination is provided in panel 9: This panel plots the end-tidal partial pressure for  $O_2$  ( $PETO_2$ ) and  $CO_2$  ( $PETCO_2$ ), and if they are measured, the exercise blood gases ( $PaO_2$  and  $PaCO_2$ ). The end-tidal curves usually progress similarly (though reversed in case of  $PETCO_2$ ) to the breathing equivalents displayed in panel 6: An initial decrease of  $PETO_2$  is followed by an upwards deflection at  $VT1$  and  $VT2$ .  $PETCO_2$  initially increases and flattens into a plateau at  $VT1$  before decreasing at  $VT2$  (and once again at  $VT3$ , provided that  $VT3$  is reached). If blood gas values are taken, the end-tidal - arterial oxygen difference  $P(ET-a)O_2$  can indicate a diffusion impairment.

Positive  $P(a-ET)CO_2$  values (which are normally negative during exercise) imply an increased dead space ventilation.

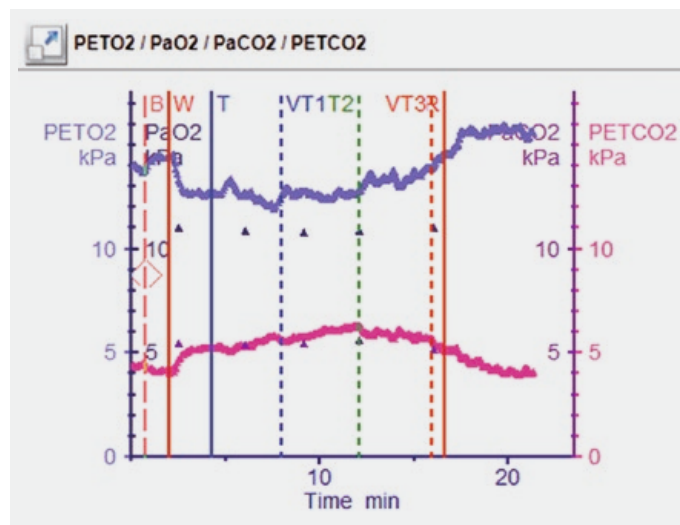

Panel 9 of the 9-Panel-Graphics.


**Attention:**

1. The  $PETO_2$  progress is similar to  $PAO_2$ . However, the alveolar gas formula is required to accurately determine the  $P(A-a)O_2$  gradient.
2. Some may approximate dead space ventilation and dead space ratio using  $PETCO_2$  instead of using blood gas values. However, this method requires caution – it may work reasonably well in healthy subjects, but in patients with certain illness, it often can provide incorrect values (Lewis (1994)).

**Energy production, Metabolism**

**Panel 8**

Conclusions regarding metabolism can be drawn based on the RER (**R**espiratory **E**xchange **R**atio, formerly known as **R**espiratory **Q**uotient RQ) by means of the panel 8, which is the ratio of  $V'CO_2$  to  $V'O_2$ . A mixed substrate metabolism of approximately 50% fat and 50% carbohydrates results in an RER value of around 0.85. A value below suggests more fat oxidation, a value above suggests more carbohydrate oxidation. In the past, the “anaerobic threshold” was usually determined by means of  $RER=1$ . Today, however, this method is no longer used: VT1 is barely recognizable by means of the RER; VT2 can usually be found near  $RER=1$ . However, this value can be only used for a rough assessment of an “upper limit” of VT2 (please see “Threshold determination” section for more details).

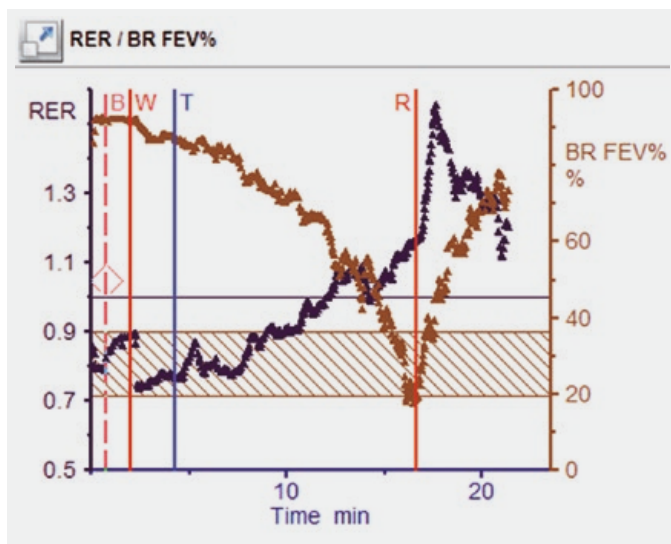

Panel 8 of the 9-Panel-Graphics.

**Note:**

In the past, RER (often also abbreviated as R) was usually termed RQ. The RQ, however, refers to the metabolism of the cell itself. Thus, the RER – which is measured at the mouth – replicates RQ only in Steady State due to the phase shift between  $V'O_2$  and  $V'CO_2$ .

**References:**

- Beaver W.L., Wasserman K., Whipp B.J.: A new method for detecting the anaerobic threshold by gas exchange. *J Appl Physiology* 60 (1986); 2020–2027
- Cooper C.B., Storer T.W.: Exercise testing and interpretation. Cambridge University Press (2001). ISBN: 0-521-64842-4
- Jones N.L.: Clinical Exercise Testing. 4th Edition, W.B. Saunders Company (1997). ISBN: 0-7216-6511-x
- Kroidl R.F., Schwarz S., Lehnigk B., Fritsch J.: Kursbuch Spiroergometrie – Technik und Befundung verständlich gemacht. 3. Auflage Thieme Verlag (2014). ISBN: 978-3-13-143443-2
- Lewis D.A., Sietsema K.E., Casaburi R., Sue D.Y.: Inaccuracy of Noninvasive Estimates of  $VD/VT$  in Clinical Exercise Testing. *Chest* 106 (1994); 1476–1480
- Roca J., Whipp B.J.: Clinical Exercise Testing. ERS Monograph 6\_2 (1997); 1-164. ISSN 1025-448x
- Rühle K.-H.: Praxisleitfaden der Spiroergometrie. 2. überarbeitete und erweiterte Auflage. Kohlhammerverlag (2008). ISBN: 978-3-17-018053-6
- Ward S., Palange P.: Clinical Exercise Testing. ERS Monograph 40 (2007). ISBN: 978-1-904097-80-8
- Wasserman K., Hansen J.E., Sue D.Y., Casaburi R., Whipp B.J.: Principles of Exercise Testing and Interpretation 3rd edition (1999) Lippincott Williams & Wilkins. ISBN: 0-683-30646-4
- Wasserman K., Hansen J.E., Sue D.Y., Stringer W.W., Sietsema K.E., Sun X-G., Whipp B.J.: Principles of Exercise Testing and Interpretation 5th edition (2012). Lippincott Williams & Wilkins. ISBN-13: 978-1-60913-899-8
- Weisman I.M., Zeballos R.J.: An Integrative Approach to the Interpretation of cardiopulmonary Exercise Testing. *Clinical Exercise Testing Prog Respir Res. Basel, Karger* 32 (2002); 300–322

# Vyntus CPX - the Latest Product Generation

## The Vyntus™ CPX at a Glance



### Vyntus™ CPX - Powered by SentrySuite™

The Vyntus CPX represents the new generation of Cardiopulmonary Exercise Testing and combines high measurement quality with ease-of-use and a workflow driven CPET evaluation. The Vyntus CPX is the result of over 50 years of experience in the development of CPET systems. The highly flexible system is suitable for various applications and can be easily used on a variety of subjects: from sick patients to top athletes and from children through to adults to old age. Thus, the Vyntus CPX can be used in a wide range of application fields. Furthermore, it ensures high-precision test results based on proven high-end sensor technology while its advanced functions offer useful support for interpretation of test results. The device is based on advanced technology and is the result of twelve generations of JAEGER, SensorMedics, ... Viasys ... Vyaire devices. It combines proven techniques with technical innovations and new medical results by offering additional measurement and evaluation procedures.



### Variable Configurations

Variable Configurations are available such as: mobile cart configuration; table top configuration; and single or dual monitor setup, and as a result the system is easily customized to your individual needs.

Combining it with a notebook will turn the system into a compact CPET station, reducing the footprint to a minimum.



**Mobile Cart configuration with dual monitor setup**



**Table top configuration with Notebook**

### Vyntus CPX covers all essential CPET applications

- Breath-by-breath Cardiopulmonary Exercise Testing
- Slow and forced spirometry, MVV as well as Pre- / Post measurements and an animated incentive
- Flow / Volume loops during exercise (EFVL) with superimposed maximum flow volume loop
- New and original 9-panel-Wasserman-graph are available along with "Possible Limitation" graph
- Ventilatory threshold determination (VT1, VT2 and VT3)
- Automatic slope calculation such as  $\dot{V}O_2$ (Watt),  $\dot{V}E(\dot{V}CO_2)$ ,  $\dot{V}O_2(\text{Log } \dot{V}E)$  (Oxygen Uptake Efficiency Slope OUES) and others.
- Ability to edit all measurement ranges for baseline, warm-up, peak, and recovery phases
- Indirect Calorimetry (REE, Fat...) using mask or mouthpiece
- Data input for RPE, blood pressure and comments
- Offline data input of blood gases with an automatic calculation of further parameters ( $P(A-a)O_2$ ,  $VD / VT, \dots$ )
- Customizable workflow for CPET evaluation
- Comprehensive program for creating individual comments and interpretations including a helpful template manager

- Automatic control of bike, treadmill and blood pressure measurement
- “CPET Protocol Editor” program for a flexible creation of ramp, step and weight dependent load protocols
- “Report Designer” program for a creation of customized reports including the possibility to export to Excel™ format at the touch of button
- “Layout Editor” to customize the graphical display during the measurement and after in evaluation mode.

## Combine Vyntus CPX with other devices or options:

- Integrated Nonin™ SpO<sub>2</sub> measurement with various sensor probes for the finger, forehead and ear
- Vyntus™ ECG: the fully integrated and wireless 12-Lead Bluetooth™ PC-ECG
- Polar™ Bluetooth™ Interface
- Choice of bike ergometers with / without integrated blood pressure measurements and treadmills
- Tango™ automated blood pressure monitor
- Indirect calorimetry using the dilution canopy method
- Measurement with high / low oxygen breathing
- Compatible with a large number of 3rd party ECGs
- Blood gas analyzer interface for serial import of blood gas data

## Optional Workflow Applications

- Questionnaire Designer and patient questionnaire application for tablets
- Networking with further PFT systems and workstations for evaluation, interpretation and central data storage
- Web-based evaluation and interpretation of PDF reports via SentrySuite Mobile Review
- Interface with hospital and medical practice systems
- Electronic Patient Records (EPR) interface through SentryConnect Interface



## Digital Volume Transducer (DVT)

The proven technology of the Digital Volume Transducer (DVT) meets the ATS / ERS guidelines for spirometry as well as passing all 24 wave forms. It is an accurate and reliable sensor for the complete flow range from low flows to maximum voluntary ventilation. Thanks to its compact and lightweight design (45 g only), the sensor has a very small dead space of only 30 mL. The DVT is insensitive both to water vapor and expired gas mix. As compared to a turbine, the flat vane system has no lag due to its small inertia. Patients and athletes will appreciate the fact that it adds minimal resistance to airflow and it is extremely comfortable to wear with both mask and mouthpiece. Different mask sizes and types (adult and pediatric, reusable or disposable) ensure best fit for each subject and ensure you can provide the highest level of clinical hygiene.




DVT with mask



DVT with mouthpiece

# Vyntus CPX

## The heart of the system - the high-precision and proven O<sub>2</sub>/CO<sub>2</sub> Analyzer

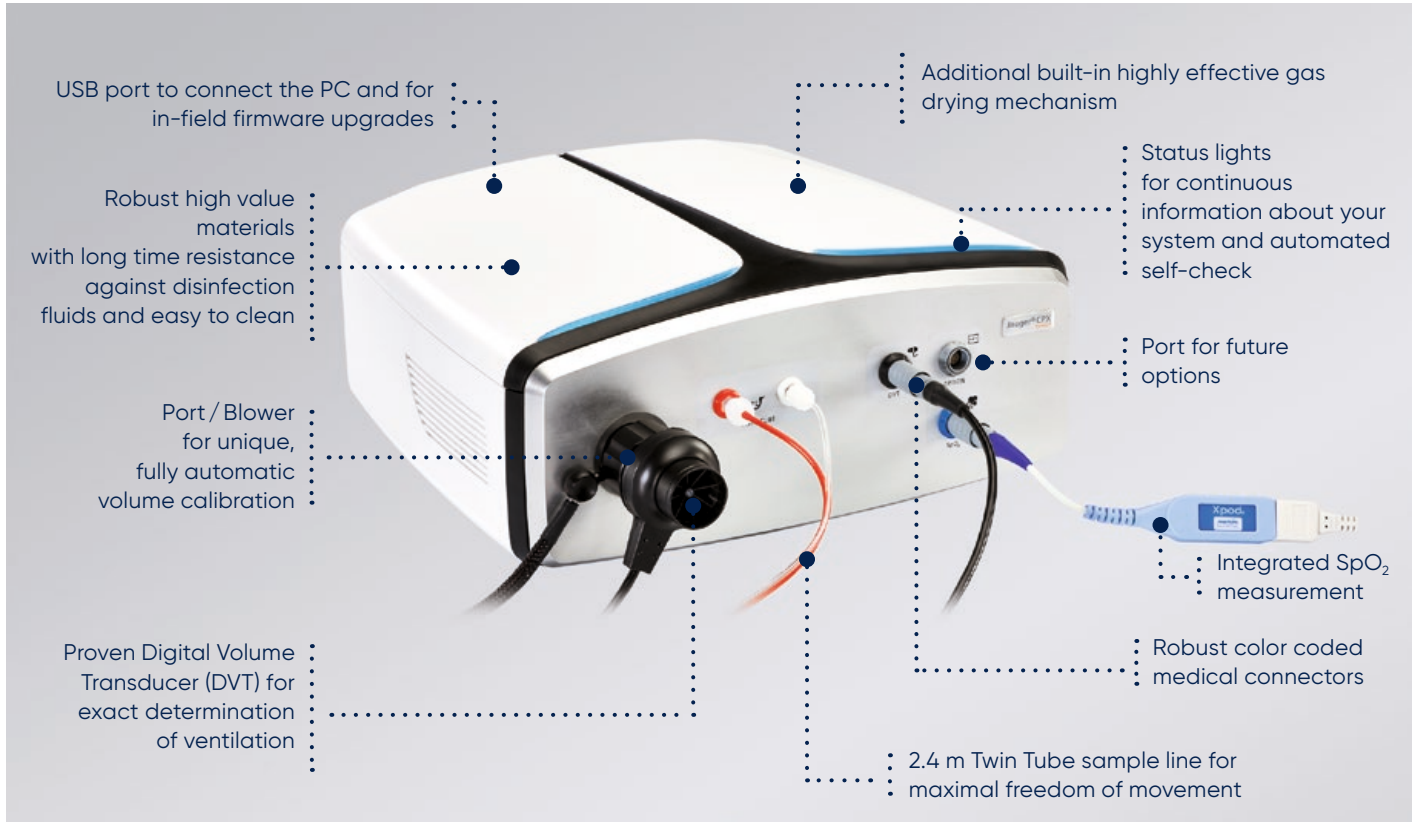

### O<sub>2</sub> cell change - made easy

The long-life (approximately 2 years) O<sub>2</sub> fuel cell can easily and quickly be exchanged at customer's site in only about a minute. All that is required is a coin to open the fuel cell door on the back of the Vyntus CPX. Take the old cell out and put the new one in.

A fully automatic filter optimization system ensures measurement continuity after the cell is exchanged.



### Calibration couldn't be easier

The Vyntus CPX is equipped with a unique, fully automatic volume calibration unit - making a manual 3 Liter calibration syringe unnecessary. Just one click in the SentrySuite

software and a volume sensor calibration will be automatically performed using the integrated blower.

The special Twin Tube (TT) sample line and the fresh air flush system allow to perform a gas analyzer calibration without disconnecting the sample tube. Additionally, the easy and fully automatic "click-and-play" 2-point gas calibration of the O<sub>2</sub>/CO<sub>2</sub> analyzers determines the delay and response times for the exact synchronization with the volume signal in one procedure.

### High accuracy and stability

Accurate and stable measurements, even during long exercise measurements, are guaranteed by: the special drying system with pre-drying via the Twin Tube; an additional arrangement inside the Vyntus CPX to remove the remaining humidity; and fast response gas analyzers (typical T<sub>10</sub>-T<sub>90</sub> = 75 msec).

### Flexibility

The Twin Tube sample line with a length of 2.4 m offers maximum freedom of movement for the patient - even with measurements performed on a treadmill.

## Vyntus™ CPX – the Software

### The Main Screen – a 360° working interface

The SentrySuite CPET Software is designed for simplified CPET testing, which can be perfectly customized to the individual needs and capabilities of your current patient. The convenient and user-friendly software interface allows easy and effective control of the measurement procedure by providing a clear overview of test and equipment controls. Furthermore, it offers valuable support for an effective interpretation of the test results. All important programs such as patient data, calibration, measurements, and even reports can be selected directly from the same screen.

The screenshot shows the main interface of the SentrySuite CPET software. At the top, a 'Patient' panel displays 'Gruber Franz' with ID 'G200367' and visit date '09.10.2018'. Callouts indicate that this panel allows users to 'Bring up personal data of the current patient, enter new patient data or search for specific patient data already available in the database'. Below this, a 'CPET BxB' dropdown menu is shown with a callout: 'Switch to other measurement programs'. To the right, a row of buttons labeled 'WS\_CPET\_E\_SHORT' is highlighted with a callout: 'Show, print and save reports'. The main window title is 'CPET BxB - SentrySuite Internal Test Version: V2.15 (SP 2.11) BETA1'. The interface includes a 'Results' table with columns for 'TIME', 'LSDC', 'HR', 'VO2RQ', 'VO2', 'VCO2', 'RER', 'VE', 'BF', 'DR FEV%\_1', and 'O2Sat'. Below the table are two graphs: 'VO2 Load / HR / VCO2' and 'DR FEV%\_1 SpO2 / Phys / HR (d)'. A 'Calibration' button with a warning icon is also visible, with a callout: 'Bring up various calibration programs'. At the bottom, a message reads 'Press F1 to start a measurement'.

Main screen after selecting the CPET measurement program. With just one click it is possible to directly switch to various programs such as patient data, calibration, spirometry, report, or to start a new measurement without leaving this window.

## The "StartUp-Window"

The screenshot shows the 'StartUp-Window' of the Vyntus CPX software. It is divided into several sections:

- Device Selection:** Shows 'Main Device' as 'Vyntus CPX' and 'Ergometer' as 'Bike ER 900'. It also shows 'Heartrate', 'Blood Pressure', and 'O<sub>2</sub>-Saturation' all set to 'Vyntus CPX'. A 'Refresh' button is present.
- Max. Predicteds:** Displays 'VE' as 123 L/min, 'HR' as 186 1/min, 'VO<sub>2</sub>' as 3074 mL/min, and 'Maximum load' as 239 W.
- PFT values:** Shows 'Measured PFT Data' with 'FEV1' at 4.05 L and 'MVV' at 8.72 L/min. 'V'E<sub>max</sub> (estimated)' is 35.0 L/min, and 'Estimated max. Load from V'E<sub>max</sub>' is 368 W.
- Profile Selection:** Shows a table of load profiles with 'ST-20-1\_(EFVL-3, BP-2)' selected as the 'Selected Profile' with a 'Maximum load' of 240 W.

Annotations on the left side of the screenshot include:

- Adjust masks and averaging settings
- Edit and create load profiles

Text descriptions for the highlighted sections:

- Device Selection:** All connected and activated devices are checked. (Green status icon: correct connection; red: incorrect connection)
- Max. Predicteds:** Automatically and individually calculated max. predicted values for the current patient. Interpretation of the final measurement results will be based on these values.
- PFT values:** Estimated maximum minute ventilation (V'E<sub>max</sub>) and estimated maximum load calculated from values measured in a prior spirometry measurement ("Measured PFT Data") to avoid the selection of an inappropriate load profile in case of a ventilatory limitation. If no measurement was performed previously, this field remains empty.
- Profile Selection:** Selection of default load profiles. Depending on the settings, the system automatically proposes the profile which comes closest to the patient's maximum load (predicted load value) - resulting from the comparison of the predicted load value with PFT maximum load - or the load profile set as standard.

The types of devices and inputs supported by the SentrySuite CPET software are divided into the categories: "Main Device", "Ergometer", "Heartrate", "Blood Pressure", and "O<sub>2</sub>-Saturation". These devices are set in the "StartUp-Window" prior to starting a measurement. In addition, the user can select from various masks and averaging methods. A wide range of preset load profiles containing many possible combinations of ergometer type (bike or treadmill) and load protocol (ramp or step load) as well as different settings for each are available from the same window. Furthermore, it is possible to create new protocols or to edit pre-existing profiles.



## Special Functions during Cardiopulmonary Exercise Testing

It is easy to control all aspects of the test procedure during the measurement by selecting the appropriate button in the left-hand button bar. Furthermore, the user can quickly switch between various displays according to their individual preference, zoom in and out each graph, and/or manually advance to the next phase.


If not already set in the pre-settings of the respective load profile, it is also possible to start an EFVL measurement, to activate the RPE scale or mark an blood-gas sampling event at any point of time during any phase with just one single click.

**F1** Next Phase  
Advance to the next phase


**F2** Profile  
Manual load change

**F5** EFVL  
Start EFVL measurement


**F6** Marker  
Mark events



9-Panel-Graphics according to Wasserman

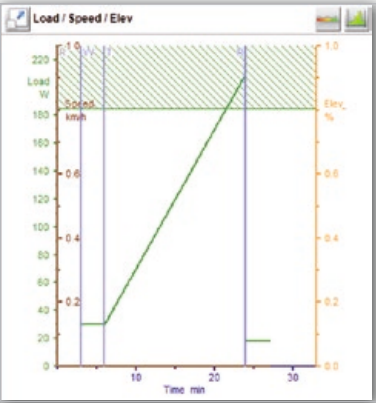

Switch from the performance graph display to the load graph display. The hatched area indicates the patient's target performance range (determined from previously calculated individual predicted values).

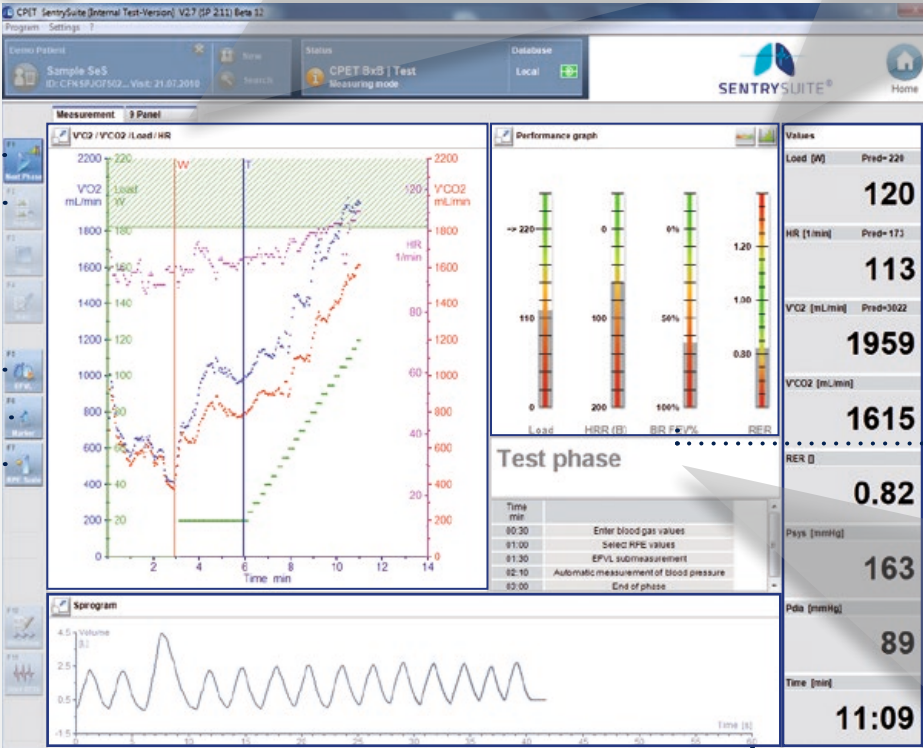

**Measurement Panel**

**Performance graph**

**Values**

Load [W]	Pred-228	<b>120</b>
HR [1/min]	Pred-173	<b>113</b>
VCO2 [mL/min]	Pred-3022	<b>1959</b>
VCO2 [mL/min]		<b>1615</b>
RER [l]		<b>0.82</b>
PaO2 [mmHg]		<b>163</b>
Pda [mmHg]		<b>89</b>
Time [min]		<b>11:09</b>

**Test phase**

Time [min]	Event
00:30	Enter blood gas values
01:00	Select RPE values
01:30	EFVL submeasurement
02:10	Automatic measurement of blood pressure
03:00	End of phase

**Spirogram**

**RPE Scale**

Exertion	Dyspnea
Extremely strong (maximal)	10
Very strong	7
Strong (heavy)	5
Moderate	3
Weak (light)	2
Very weak	1
Extremely weak	0.5
Nothing at all	0

Currently measured values (predicted and actual values) in numbers

Performance graph

**Test phase**


Time [min]	Event
01:00	Automatic measurement of blood pressure
02:30	EFVL submeasurement
03:10	Automatic measurement of blood pressure
05:10	Automatic measurement of blood pressure

This table provides a convenient overview of the measurement process. Upcoming events including their times are displayed enabling the user and the patient to get prepared. The displayed time is a countdown showing time (in minutes and seconds) remaining for the event to take place. The displayed events and phase duration depend on the settings in the selected load profile.


# Vyntus CPX Software

## Everything at a glance

SentrySuite optimizes test efficiency by means of concise graphical overview. The performance graph clearly displays the degree of the maximum load with regard to "Load", "Heart Rate Reserve (HRR (B))", "Breathing Reserve (BR FEV%) and "Respiratory Exchange Ratio (RER)":

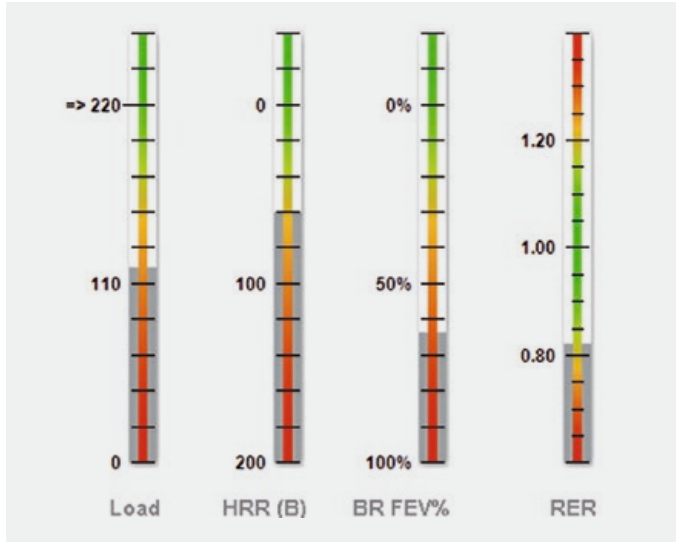


The filling bars indicate the current performance of the subject (grey areas). If large reserves are available, the load can correspondingly be increased.

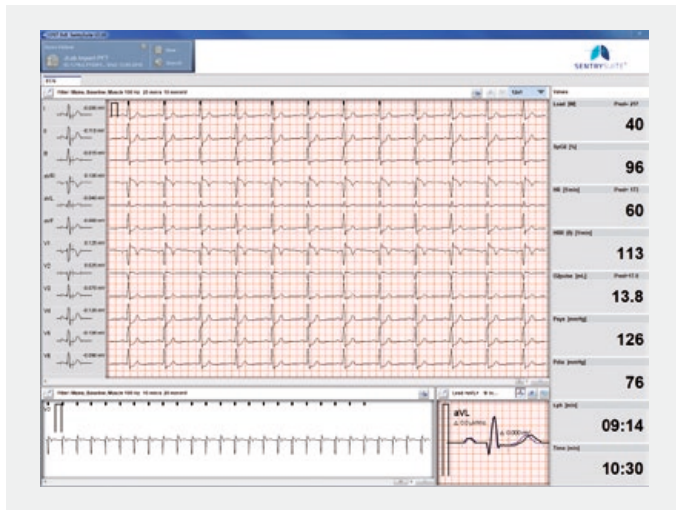

It is also possible to record heart rate (with or without ECG) during the exercise test. With its fully integrated 12-lead Vyntus™ ECG for rest and stress ECG, Vyaire Medical offers the optimum solution for this purpose: a complete ECG recording on a second monitor or as a single graph (using a single monitor). Alternatively, the heart rate can also be recorded via the integrated SpO<sub>2</sub> sensor, a Polar™ chest strap or other, combinable 3rd party ECG systems for a comprehensive CPET measurement.

The numerical display of parameters allowing the user to easily and quickly read both actual and predicted values is perfectly supplemented by the graphical display of the respective parameters.

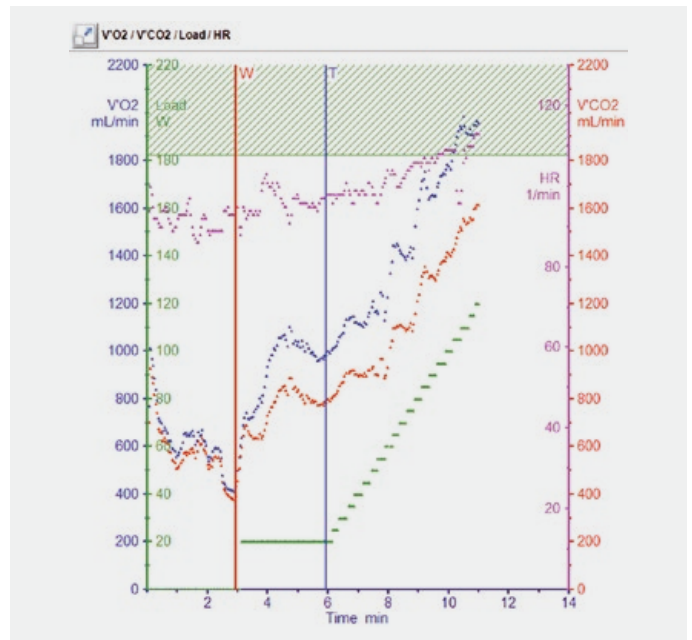

Graphical display of the actual and predicted value(s)

In this graphic, the user is able to see the parameter values V'O<sub>2</sub>, V'CO<sub>2</sub>, load and HR from the current phase as well as from the previous phases at a glance. The vertical lines indicate the start of the respective phase. Additionally, the hatched area indicates the predicted target load range for a normal subject or patient.

Value	
Load [W]	Pred=220 <b>120</b> Actual and predicted load value (Watt)
HR [1/min]	Pred=173 <b>113</b> Actual and predicted heart rate
V'O <sub>2</sub> [mL/min]	Pred=2022 <b>1959</b> Actual and predicted oxygen uptake
V'CO <sub>2</sub> [mL/min]	<b>1615</b> Actual carbon dioxide production
RER [ ]	<b>0.82</b> Actual "Respiratory Exchange Ratio"
Psys [mmHg]	<b>163</b> Actual systolic blood pressure
Pdia [mmHg]	<b>89</b> Actual diastolic blood pressure
Time [min]	<b>11:09</b> Elapsed time (total)

Numerical display of the actual and predicted values. It is possible to customize both parameter selection as well as parameter sequence.

## Results and Interpretation

The ability to switch between various graphical displays in the result screen provides the user both assistance and a selection of various approaches for interpretation. Many of the displays may also be customized to the users own preferences.

9-Panel-Graphics according to Wasserman (2012)

EFVL measurement results

BR/HRR against time

Measurement results displayed numerically

Possible Limitations

Workload graph

Measurements can be reviewed. You can use the buttons in the patient field to select patients and visits.

The "Guidance" tab provides textual assistance

$V'O_2$ ,  $V'CO_2$ , HR and Load against time

Automatic interpretation and classification of the measurement results

User comments/interpretation

Time min	Load W	HR min	VCO2g mL/min/kg	VE mL/min	VCO2 mL/min	RER	VE L/min	BP mmHg	BRFEV% %	O2uptake mL
00:14	0	80	1.6	136	120	0.82	6	6.3	94	2.3
00:24	0	81	3.5	325	277	0.82	6	6.4	94	4.9
00:37	0	80	2.9	281	224	0.89	6	6.9	97	4.2
00:47	0	80	3.9	354	298	0.85	6	6.7	94	5.1
00:52	0	80	3.7	319	287	0.84	6	6.5	96	5.0
01:08	0	81	3.8	323	277	0.87	6	6.3	97	5.6

# Vyntus CPX Software - Evaluation

## Graphical Display of the Results

In the **Result** screen, the data is displayed numerically as well as graphically, and can be shown as either time or breath averaged. It is also possible to select additional displays in the right-hand graphic (BR / HRR).

**Results**

The left window displays:  $\dot{V}O_2$ ,  $\dot{V}CO_2$ , Heart Rate (HR) and Load (Watt) against time

Time min	Load W	HR $\dot{V}$ min	$\dot{V}O_2$ /kg (mL/min)/kg	$\dot{V}O_2$ mL/min	$\dot{V}CO_2$ mL/min	RER	VE L/min	BF 1/min	BR FEV% %	O2pulse mL
00:14	0	60	1.6	136	120	0.88	5	6.0	98	2.3
00:24	0	61	3.5	302	277	0.92	9	5.6	96	4.9
00:37	0	60	2.9	251	224	0.89	8	5.9	97	4.2
00:47	0	60	3.5	304	258	0.85	9	5.7	96	5.1
00:52	0	60	3.7	318	267	0.84	9	5.5	96	5.3
01:06	0	61	3.6	309	236	0.77	8	6.3	97	5.1

The right window displays: Breathing Reserve (BR FEV%) and Heart Rate Reserve (HRR) against time

The hatched areas mark the individual predicted areas of the subject

The right window displays: Breathing Reserve (BR FEV%) and Heart Rate Reserve (HRR) against time

If you move the dashed line (with diamond in the centre) to the left or to the right in one of the two windows, the values displayed in the tabular data will change according to the new position. If combined with Vyntus ECG, ECG tracing will also move to newly displayed position.

The vertical lines in the charts display the different markers during the measurement; e.g.:

- W** = Start of the warmup phase
- T** = Start of the exercise phase
- R** = Start of the recovery phase

It is possible to add additional markers.

The points VT1, VT2 and VT3 define the ventilatory thresholds. (For more information over the ventilatory thresholds please see the "Threshold Determination" chapter in this special edition.)

# Vyntus CPX Software - Evaluation




Besides the aforementioned display of predicted values, numerical and graphical measurement values, and 9-Panel-Graphics, the software offers further valuable graphics and tools which significantly simplify and support the evaluation of results. These additional features are:

## Possible Limitations

The possible limitation chart indicates and excludes possible and specific diseases as the cause of cardiopulmonary limitations. However, it must be kept in mind that the possible limitations are only a suggestion based on the measurement data and need to be verified by the user (modified according to Weisman (2002)). The SentrySuite CPET Software checks the measurement results for the following possible limitations:


- Heart failure
- COPD
- Interstitial lung disease
- Pulmonary Vascular disease
- Obesity
- Deconditioning

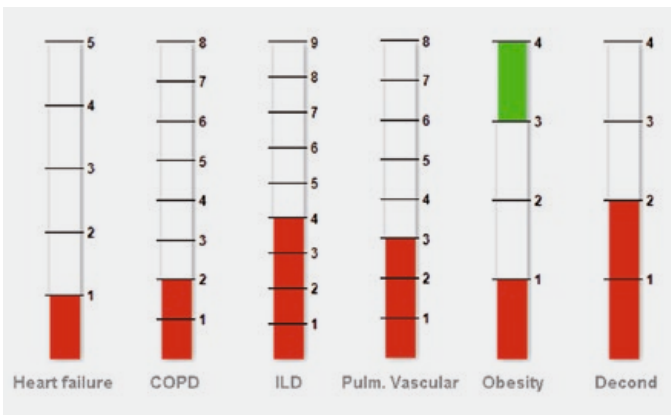


Every visible segment in a bar corresponds to a parameter. If the parameter reaches a limit value which is indicative of a limitation, the respective segment will be highlighted in red. If a limitation can be excluded due to the measured values, the respective segment in the respective limitation bar will be marked green. White segments indicate normal values whereas potential grey highlighted segments point out that an evaluation of the respective parameter is not possible, usually due to missing blood gas values (for example to calculate  $P(A-a)O_2$ ).

## EFVL - Flow/Volume during exercise

Subjects with limited lung function can only be subjected to physical exercise to a limited extent and therefore need to be observed carefully during a CPET measurement. An EFVL measurement (Exercise Flow Volume Loop) allows to supervise these patients and to decide whether to continue the exercise or terminate the test. It also indicates important aspects on when or whether a lung function disorder prevents further exercise.

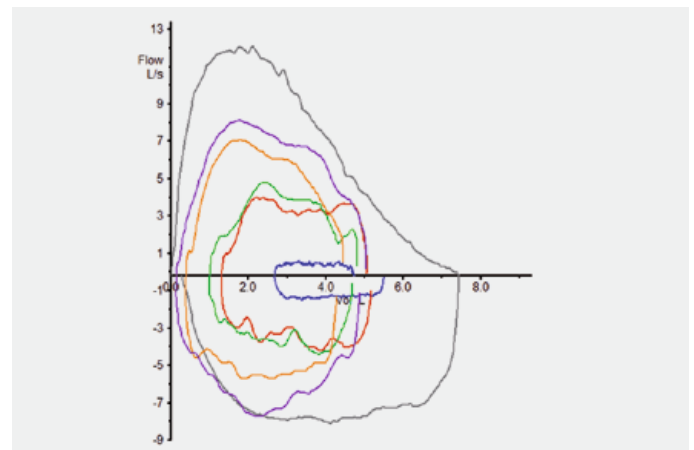

Measured Flow-Volume loops (superimposed)

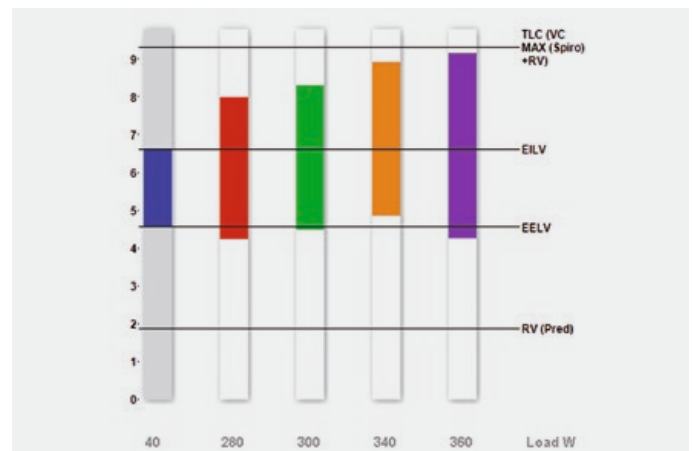

Flow-Volume Curves recorded during CPET displayed as a bar diagram

# Vyntus CPX Software - Evaluation


## Vyntus CPX Evaluation Workflow - easy to use from beginners to experts

After the measurement, the evaluation workflow will automatically guide you step-by-step through post-test editing (also referred to "EditMode"). Just click "Next" to move from one step to another. This procedure standardizes your evaluation / interpretation and reduces your time to produce a result. Depending on the use, it is possible to create different workflows for different users or user groups. The program will support you with threshold determination or slope calculation by means of intelligent evaluation. The final decision however, is up to the interpreting clinician.

The entire workflow includes

- Entry of End of Test Criteria, manually or by means of predefined texts
- Editing the ranges of rest, warm-up, test and recovery phase
- Editing the ranges for slope determination purposes
- Editing the ventilatory threshold VT1
- Editing the ventilatory threshold VT2
- Editing the ventilatory threshold VT3
- Editing the measured EFVL (Exercise Flow / Volume Loops)
- Editing RPE / Entering or editing markers, blood gases, RPE values ...
- Editing Steady State measurements

### Example Workflow:

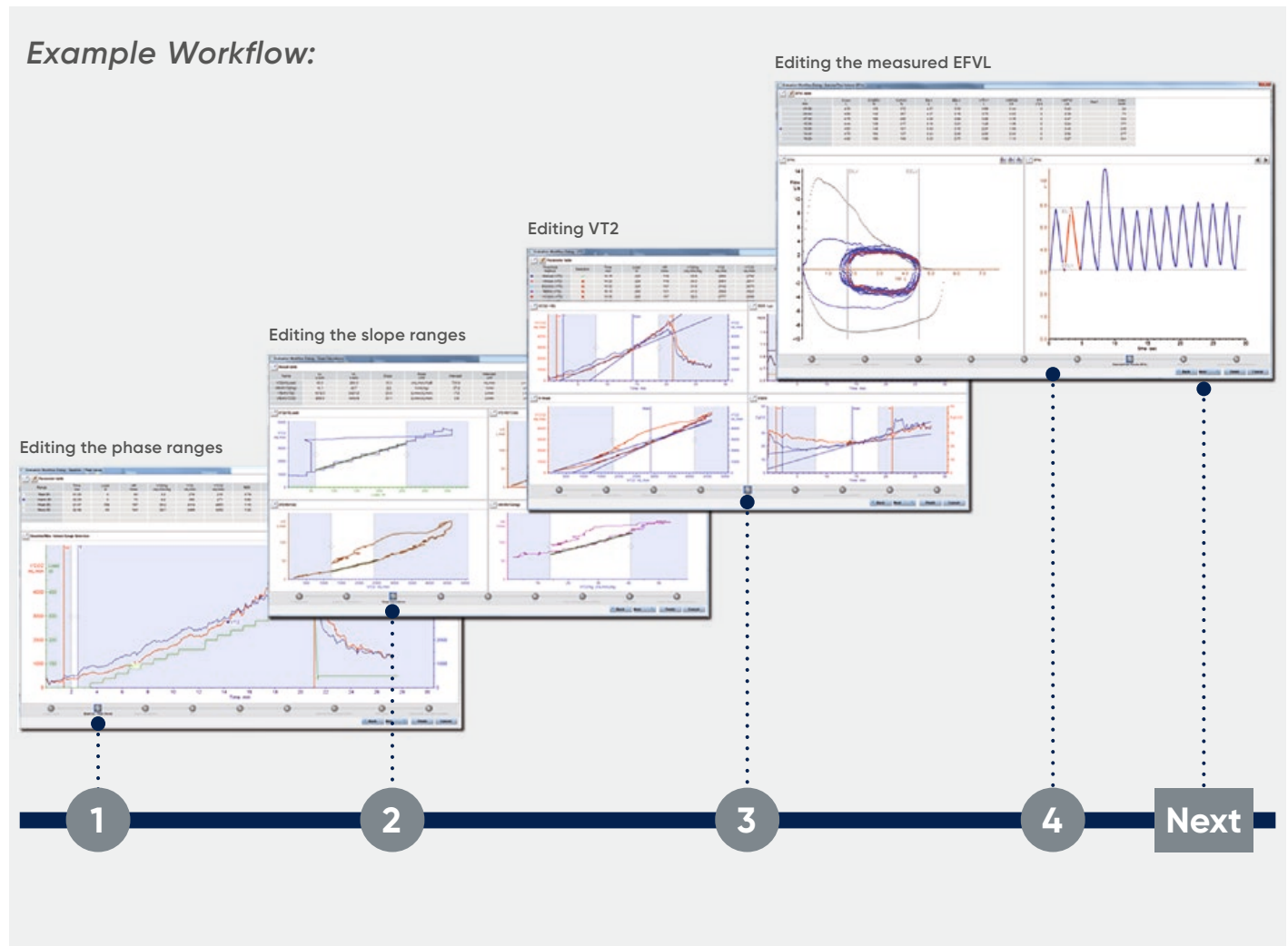

## Automatic Interpretation and Classification of the Measurement Results

The "Auto Interpretation" tab provides an automatic textual interpretation of the measurement results:

In addition to the textual interpretation, a classification of the test results is displayed. The classification is based on the predicted value of maximum oxygen uptake (Löllgen (2010)):

Excellent	=	$V'O_2\text{max \% Pred}$	$\geq$	120
Normal	=	$85 \leq V'O_2\text{max \% Pred}$	$<$	120
Mild	=	$70 \leq V'O_2\text{max \% Pred}$	$<$	85
Moderate	=	$50 \leq V'O_2\text{max \% Pred}$	$<$	70
Severe	=	$V'O_2\text{max \% Pred}$	$<$	50

... With Auto Interpretation, it is possible to choose between several authors for a suggested interpretation. The respective measurement program saves the selected author as the standard author for the next examination. Among others, the authors "CPET Eschenbacher, Mannina (1990)" - Eschenbacher (1990) are available. Note that auto-interpretation does not substitute for medical advice, provides only support for qualified personnel and shall always be reviewed by a physician.

## User Comments/ Interpretation

The "Interpretation/ Comments" tab allows the user to enter individual comments and/or interpretations manually. It is possible to load various templates and macros or to customize texts entered manually by choosing between various layout features. This allows the flexibility to quickly create comprehensive customized reports on the day of the test. This can be taken to the extent where patient clinic letters can be completed within the software without having to resort to additional dictation and letter writing.

**Technician**  
 A cardiopulmonary stress test was performed on **Demo CPET-Jose** using Carefusion Vyntus CPX metabolic system. The patient exercised on a cycle ergometer using a ramp protocol. Heart rate was measured using the Vyntus ECG system.  
 The results of the test are as follows:

**Physician**  
**FUNCTIONAL CAPACITY:**  
 The patient exercised for **9** minutes. The peak work rate achieved was **196 W**. The maximum oxygen uptake (VO2 Peak) was **2263 mL/min, 90 %** of predicted. The VO2max/kg was **25.1 (mL/min)/kg, 89.8 %** of predicted.  
 The ventilatory threshold (VT1) occurred at **1114 mL/min, 44 %** of predicted max VO2 (normal >40%).  
 The ventilatory threshold (VT2) occurred at **1590 mL/min, 63 %** of predicted max VO2 (normal > 60%).

Choose standard text modules under "Templates" or compose an individual text. According to the template, the corresponding measurement values will be imported directly from the measurement and incorporated into the text.

Thus, an entire summary can be created with one single click, which - if necessary - can be edited or extended. Both graphics as well as measurement and evaluation parameters can be transferred to predefined or user-generated reports. Easy export of the data into Excel™ for further processing is also possible.


# Vyntus CPX Software - Edit Mode

## Edit Mode

It is not only possible to edit the automatic workflow at the end of a measurement, but single sections may also be edited by means of the "Edit Mode", which shall only be described as an example in the following.

### Threshold determination

As pointed out in the later section "Threshold Determination" (see chapter "Basics and Diagnostics"), there are different procedures for the determination of the respective thresholds. These shall be discussed by means of VT1:

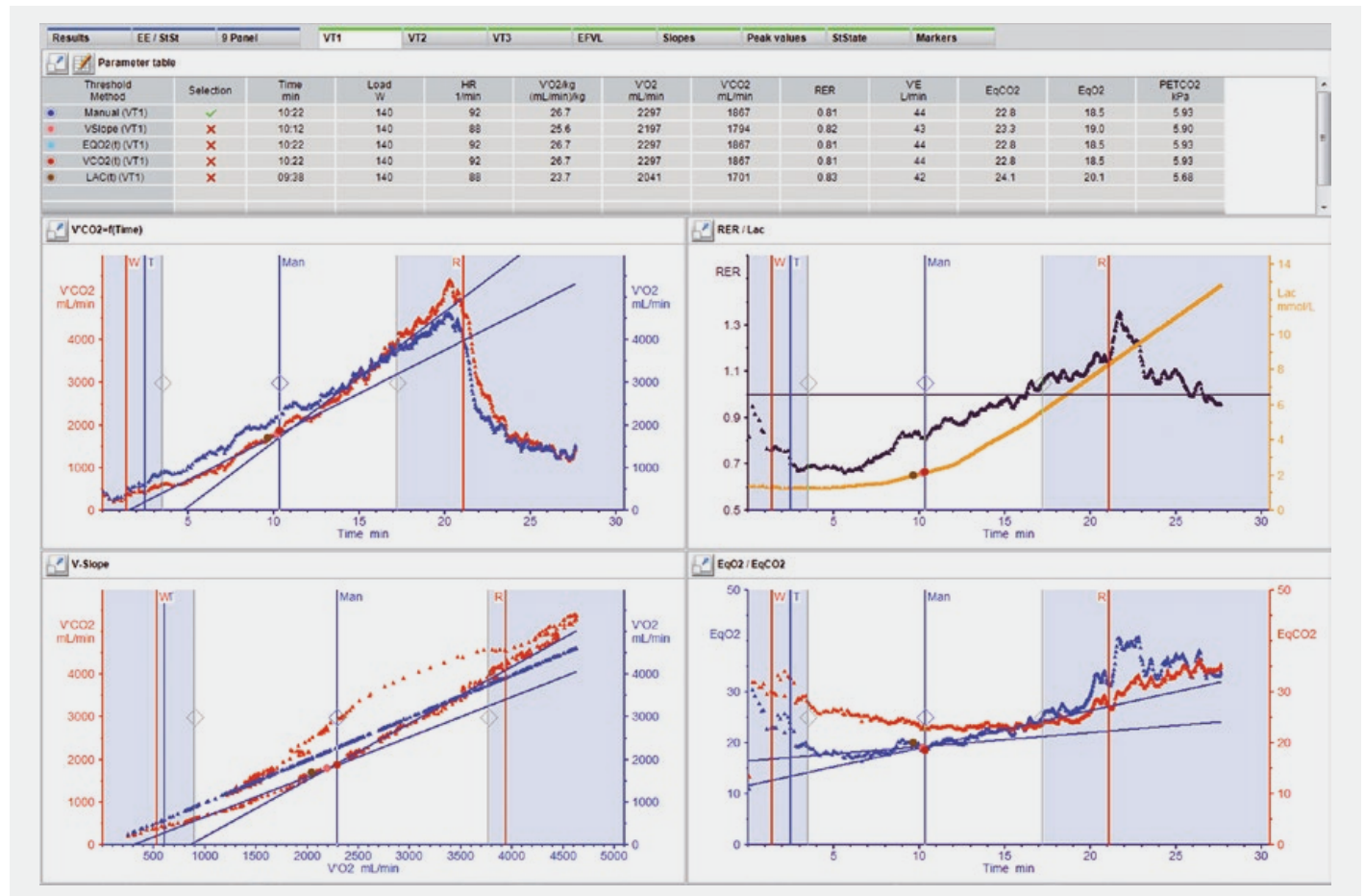

At the end of a measurement the program tries to mathematically determine the different thresholds (break points) within the specified white area and marks them:

- Orange: Break point in the V-Slope Graphic ( $V'CO_2(V'O_2)$ )
- Light blue: Break point in  $EqO_2(\text{Time})$  - in this example superimposed by the red point.
- Red: Break point in  $V'CO_2(\text{Time})$

As those break points are usually not identical, the average of all determined break points is displayed (vertical blue line). Furthermore, the program tries to confirm the break points by means of the regression line.


In case the users do not agree, they are free to modify the white areas to initiate a recalculation or alternatively manually shift the blue line to the position they believe the threshold to be. The corresponding data will also be displayed numerically in the table at the top. For a better evaluation, each graph can be expanded to full screen just by one mouse click.



## Slope Calculation

As already outlined in the introduction, various evaluations also require the calculation of the dynamic behavior of parameters (e.g.  $V'E(V'CO_2)$ -Slope).

This calculation is automatically performed already at the end of a measurement. Via the edit mode, it can be checked, and if applicable, edited.

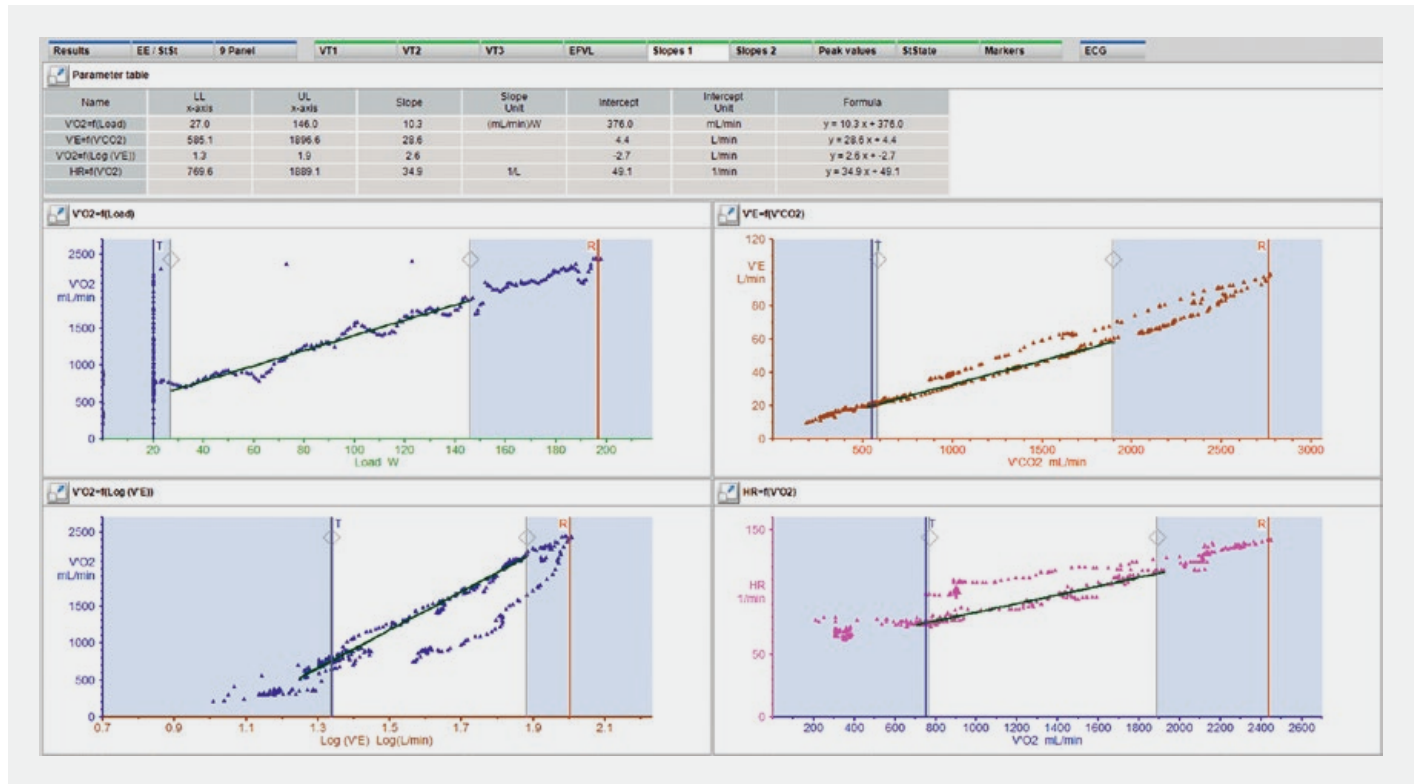

The 4 most important slopes calculated by means of the white areas are the following:

- Top left: Aerobic capacity ( $V'O_2$ (Watt))
- Top right: Respiratory efficiency for  $CO_2$  ( $V'E(V'CO_2)$ )
- Lower left: Respiratory efficiency for  $O_2$  (OUES:  $V'O_2$ (Log(V'E)))
- Lower right: Cardiovascular efficiency ( $HR(V'O_2)$ )

In this mode, the user can modify the pre-set white areas and thus initiate a recalculation.

The corresponding data will also be displayed numerically in the table.

### References:

- Eschenbacher W.L., Mannina A.: An algorithm for the interpretation of cardiopulmonary exercise tests. Chest 97 (1990); 263-267.
- Löllgen H., Erdmann E., Gitt A.K.: Ergometrie, 3. Edition. SPRINGER (2010).
- Weisman I.M., Zeballos R.J.: An Integrative Approach to the Interpretation of cardiopulmonary Exercise Testing Clinical Exercise Testing Prog Respir Res. Basel, Karger 32 (2002); 300-322.

# Vyntus CPX - Options

## Vyntus™ ECG

The Vyntus ECG is wirelessly connected to a PC or notebook via Bluetooth™ and allows for a recording of a 12-lead resting as well as stress ECG.



### Vyntus ECG - the art of diagnostic integration

The Vyntus ECG is intended for measuring the surface ECG of the patient. It communicates wirelessly and directly via Bluetooth™ and integrates conveniently with the Vyntus CPX system. Patients will appreciate the wireless technology, the small and light design of the amplifier and the short electrode cables which improves comfort and providing maximum freedom of movement. If an untoward event occurs during testing, the wireless connection of the ECG permits easier movement of the patient to a table or chair, while maintaining constant ECG collection and display. Additionally, the all-in-one view ensures a user-friendly interface. The acquired ECG can be displayed on the screen or conveniently be printed on paper.

### Excellence in diagnostic and prognostic value in a powerful combination

Exposing the heart to increased workloads is often the only way to detect cardiac abnormalities. Consequently, the combination of heart and lung parameters is essential for comprehensive cardiopulmonary exercise testing. The Vyntus ECG allows for a 12-lead stress ECG recording while automatically evaluating and analyzing the signals. Detected abnormalities such as extrasystoles or pacemaker control are displayed on screen during the measurement. The user can modify

- Speed
- Gain
- Lead selection
- Define print areas
- Print an on-line report during the measurement

Additionally, the Vyntus ECG provides a "Full Disclosure" feature for saving the unfiltered, continuous ECG signals. Via SentrySuite, the 12-lead Vyntus ECG integrates fully and seamlessly with the Vyntus CPX system. This enables laboratories to leverage their medical devices as well as Healthcare IT investments and provides an easy and clear interpretation of the measurement results.


Further benefits:

- One user interface
- One program to train
- One central database
- One combined report
- One network interface
- One HIS-connection



## ECG Recording

As soon as all electrodes are connected, the minimal fast potential differences originating from the heart can be detected on the body surface and subsequently be recorded by the Vyntus ECG. At the beginning of the measurement, the electrode contacts are checked automatically.

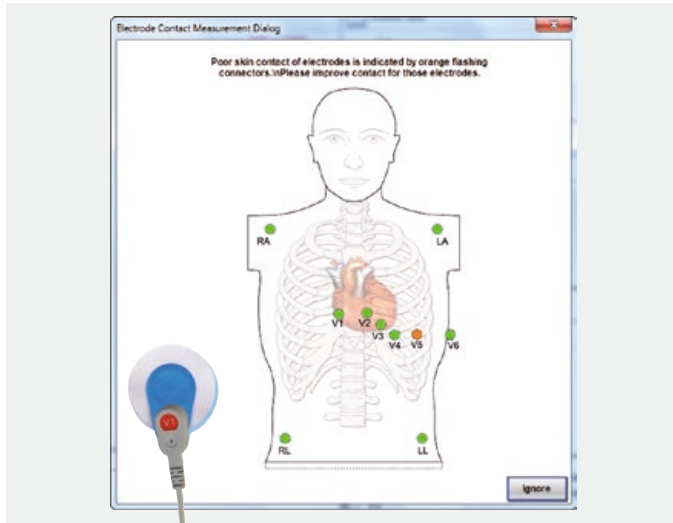

The green electrodes displayed in the screen indicate good contact. If a contact is poor, the respective electrode flashes orange.

The ECG recording starts automatically as soon as the electrodes and the skin are in good contact.

## Resting ECG


If required, several resting ECG trials can be recorded and compared to each other (similar to spirometry). In addition, a proposed interpretation according to HES (Hanoverian ECG Interpretation system, Willems (1991)) is compiled.



Resting ECG after recording with an interpretation according to HES as well as a classification

## Stress ECG


The stress ECG application offers an attractive graphical user interface and leaves nothing to be desired:

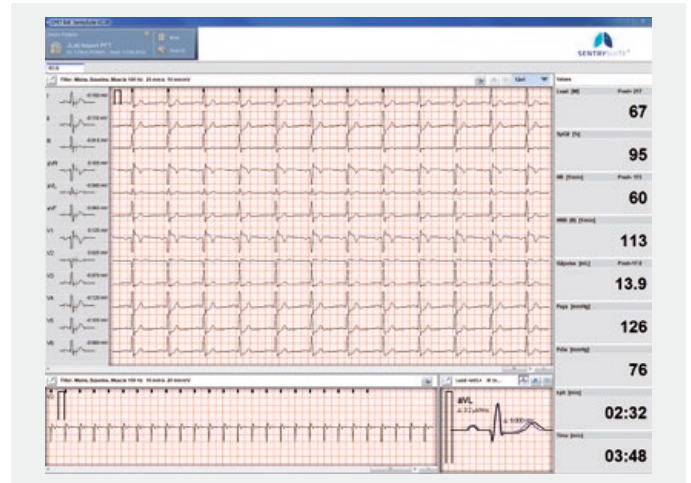

Stress ECG during the measurement

As well as the continuous recording of the single leads, the complexes, including the appropriate ST values, are displayed on the left. The lower left area of the screen shows a full disclosure recording with potentially present abnormality markers. Both recordings can be paused and scrolled back during the measurement offering a close look at previous signals. The complex shown at the lower right screen section is displayed with the reference signal including the appropriate numerical abnormalities and can be customized by the user for speed, gain and lead selection.

**All ECG raw data is recorded and saved during the entire exercise measurement.**

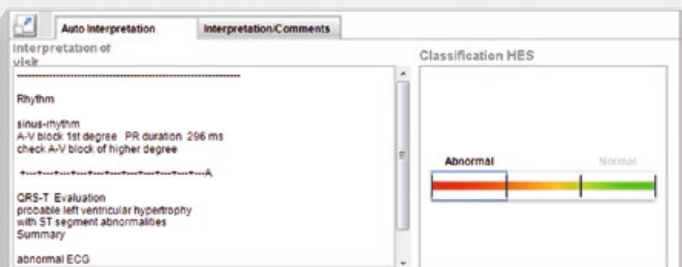

The HES™ program was part of the “Common Standards for Quantitative Electrocardiography” project, CSE. The results were independently analyzed in


**Willems J.L et al.:** The diagnostic performance of computer programs for the interpretation of electrocardiograms. N Engl J Med. 25\_325 (1991); 1767-73.

# Vyntus CPX - Options

## Vyntus™ CPX - High / Low FIO<sub>2</sub> Option (Not available in US)

### A powerful extension

This option allows the user to make measurements while the subject breathes increased or decreased concentrations of inspired oxygen. For this, a Y-valve is connected to the volume sensor (which is also connected to the gas sample tube) permitting the subject to inhale the prescribed oxygen concentration from a reservoir and allowing a CPET measurement to be performed simultaneously.

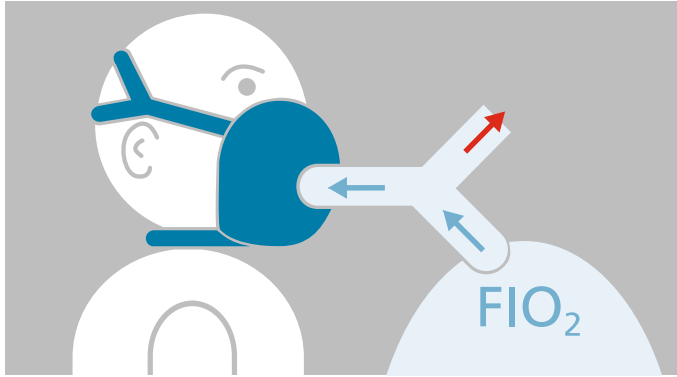

Measurement principle; the breathing bag containing the oxygen concentration to be inspired can either be refilled via a gas cylinder or by means of an appropriate blender.



Arrangement of the individual parts; the measurement can either be performed with a mask or with a mouthpiece. A head-gear for support is available as well.

### High FIO<sub>2</sub>

Subjects suffering from a ventilation-perfusion disorder (e.g. transplant patients, idiopathic pulmonary fibrosis, severe COPD) are often not able to handle everyday life without supplemental oxygen. In order to at least perform some simple tasks, such as moving around home or taking a walk, they are often equipped with a portable nasal oxygen supply. To examine those subjects' exercise capacity, the patient must be supplied with additional oxygen during the measurement. However this is not possible using nasally supplied oxygen, because as ventilation increases with exercise, it would have the effect of diluting the oxygen and consequently lowering FIO<sub>2</sub>. In this situation wash in or wash out effects would cause the user to measure the superimposition of oxygen uptake and wash out effects (or wash in effects, respectively) rather than the actual oxygen uptake. In order to avoid this problem, the subject is provided with a constant FIO<sub>2</sub> concentration (typically 30% – 40%) via the breathing bag during the measurement. The measurement procedure is otherwise similar to a normal BxB measurement. Additionally, this option considers the additional dead space caused by the Y-valve. With these measurements, special attention must generally be paid to the prior washing in of the lung and the blood (this is evident when V'O<sub>2</sub> is too high and RER is consequently too low) until balance is reached. If a severe ventilation-perfusion disorder is present, this can take up to 10 minutes. Only then the actual measurement and exercise should be started.



Behind the scenes, however, the software applies the Eschenbacher transformation (Eschenbacher (2016)) for the calculations as the Haldane transformation (Haldane, (1912)) does not provide plausible and reliable data, especially at high FIO<sub>2</sub> values (also refer to chapter "Haldane and Eschenbacher transformation").

## Low FIO<sub>2</sub>

The same procedure can be used to reduce the inhaled oxygen “concentration”. The oxygen uptake strongly depends on the oxygen partial pressure PAO<sub>2</sub> in the lungs and therefore on the environmental partial pressure PIO<sub>2</sub>.



Low FIO<sub>2</sub> used to simulate high altitudes

The current air pressure Pbar and the oxygen concentration FIO<sub>2</sub> result in

$$PIO_2 \text{ [kPa]} = Pbar \text{ [kPa]} * FIO_2 \text{ [%]} / 100$$

Consequently, the partial pressure of environmental oxygen (and, with it, the oxygen partial pressure in the lungs) strongly depends on the environmental air pressure and thus on the altitude.

From the above equation a low PIO<sub>2</sub> can be achieved by reducing either the air pressure or the inspired oxygen concentration. The relationship between altitude and FIO<sub>2</sub> can easily be estimated and is shown in following table:

Altitude [m]	Pressure [hPa]	FiO <sub>2</sub> [%]
0	1013	20.9%
500	956	19.7%
1000	901	18.6%
1500	849	17.5%
2000	799	16.5%
2500	752	15.5%
3000	707	14.6%
3500	664	13.7%
4000	624	12.9%
4500	585	12.1%
5000	549	11.3%
5500	514	10.6%
6000	481	9.9%
6500	450	9.3%
7000	420	8.7%
7500	392	8.1%
8000	366	7.6%
8800	326	6.7%

Correlation between altitude, air pressure and oxygen concentration

Instead of performing a measurement at an altitude of 2500m, an identical PIO<sub>2</sub> can be established in a laboratory by reducing the FIO<sub>2</sub> (at sea level) to approximately 15.5%. Doing this offers the ability to simulate altitude and examine a patient’s response to altitude within the safe environment of a laboratory, for example: will they desaturate on an airplane during a flight; or will an athlete respond well to the effects of high altitude training.



A modern passenger plane has a pressure complying with about 2500 m (or ca. 15.5% O<sub>2</sub>).

The measurement procedure is almost identical to a BxB measurement with an increased oxygen supply; the only difference being that the breathing bag is filled with a reduced oxygen concentration.

Again, special attention must be paid to the prior washing in of the lung and blood (V’O<sub>2</sub> too low, RER too high) until equilibrium is reached. Only then the exercise should start and accurate measurement can be taken.

Furthermore, it should be noted that with both a high and low FIO<sub>2</sub> the adjusted oxygen concentration needs to remain constant during the entire measurement in order to avoid wash in and wash out effects during the measurement.

### References:

- [1] Haldane J.S.: Methods of air analysis. Charles Griffin & Co., Ltd., JB Lippincott Co., Philadelphia (1912).
- [2] Eschenbacher H.: Haldane and Eschenbacher transformation. White Paper RC8168 (0118/PDF). Vyair (2018).

# Vyntus CPX - Options



## Vyntus™ CPX - Canopy Option (Not available in US)

### Canopy - Indirect Calorimetry

Indirect calorimetry measurements can be made with a face mask or a mouth piece using breath by breath technology, however the dilution canopy option improves subjects' comfort: This dilution canopy method is a proven and patient-friendly method to determine the resting energy expenditure. During the development of this option, increased focus was placed on our cleaning and hygiene concept, which is becoming a progressively more important topic within the hospital and laboratory environment.

By means of an innovative, patented dilution system including a single use canopy (1), the plug-in rings (3) used to attach the transparent foil to the holder (2) are the only parts which need to be cleaned and disinfected between the measurements.

Parts of the powerful Canopy module (6), used to draw the sample gas through the system, can be disassembled and cleaned (7) but the design allows for the additional application of a bacterial filter (4) to avoid contamination of all downstream parts (5, 6 and 7) including the DVT.

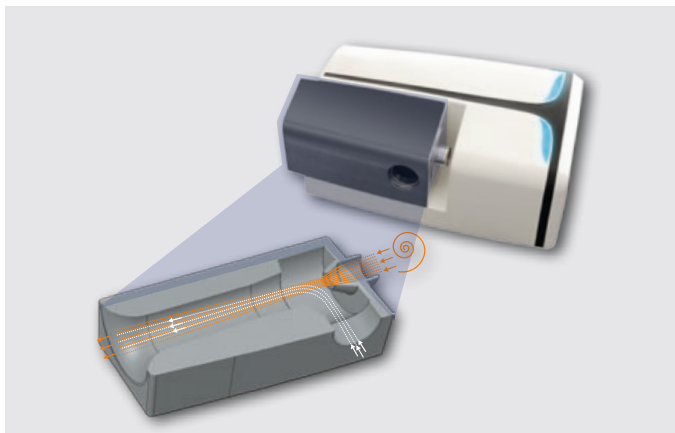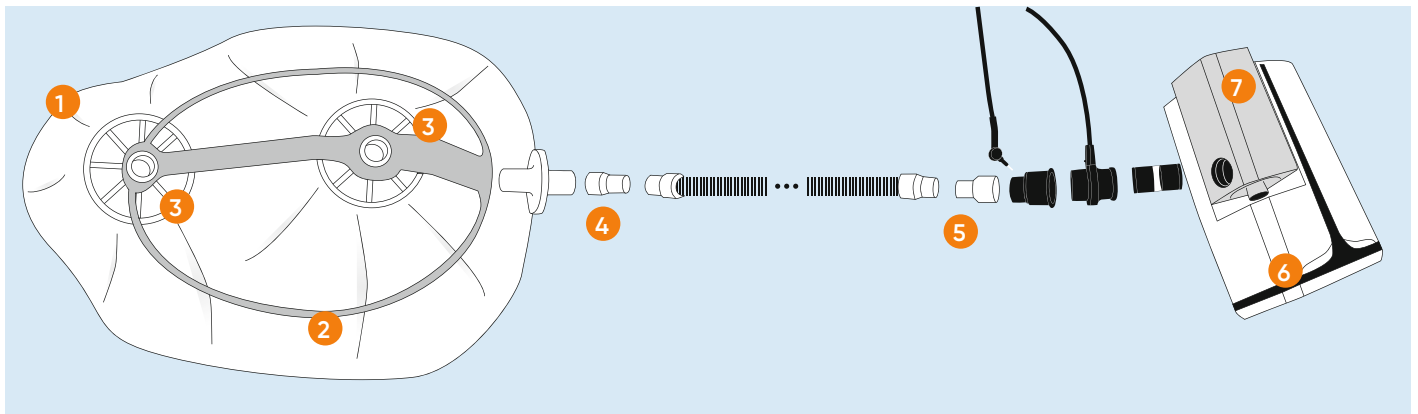


The dilution system flow can be varied over a wide range reaching from approx. 25 L/min to approx. 80 L/min.

It can also be adapted to a specified value or automatically be controlled by the software. Thus, the dilution flow can automatically be adapted to the patient by means of the measured CO<sub>2</sub> concentration.

The blower-system is constructed so, that the exhaled air from the patient does not come into contact with the blower itself.

## Measurement and Evaluation

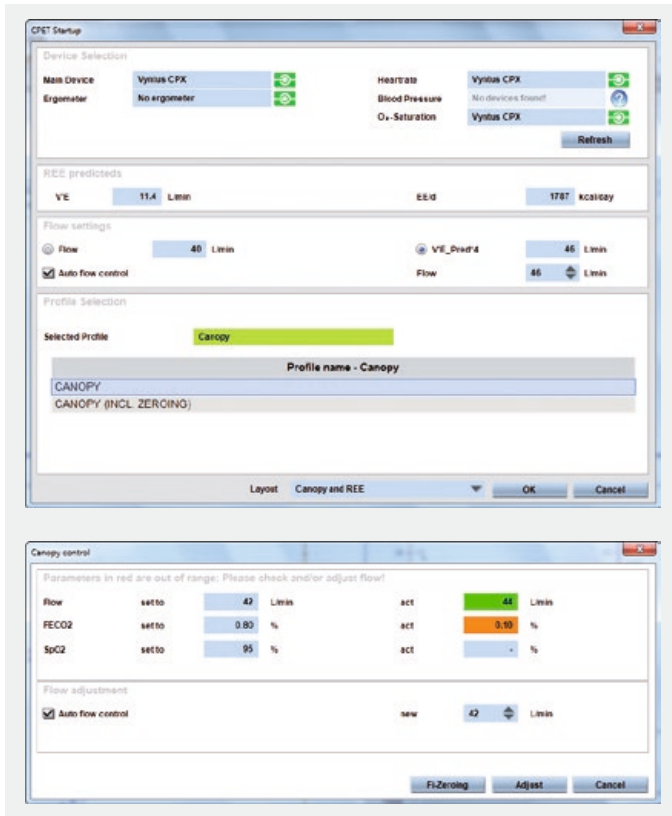

The measurement procedure for Canopy mode is similar to a BxB measurement. However, the CPET-startup window also shows the predicted values for resting energy expenditure as well as fan control.

In the section beneath, the SentrySuite software uses patient data to propose the initial flow to be used for the measurement. Additionally, the user can select whether to perform the measurement using a specified value for V'E or if the software should control the flow automatically (these settings can also be modified during the measurement).

As with BxB, you can select various profiles to automatically run additional features such as: zeroing of the gas analyzers during longer measurements; taking blood pressure; or blood gas entry during the test.

During the measurement, flow, FE<sub>CO<sub>2</sub></sub> and SpO<sub>2</sub> are constantly monitored and controlled. If a deviation from the predicted values is detected, a window will point this out and offers the user the option to intervene.

As soon as the measurement is finished, the data is clearly displayed both numerically and graphically. The user has the ability to discard or summarize single areas for evaluation purposes. If required, the measured urinal nitrogen value can be entered for further calculation.

Among the usual parameters such as V'O<sub>2</sub>, V'CO<sub>2</sub>, RER and V'O<sub>2</sub>/kg, the following parameters are also of interest:

- EE (Energy Expenditure)
- npRER (non protein RER)
- Division of energy production and substrate consumption in
  - Carbohydrates
  - Fats
  - Proteins (with the entry of urinary nitrogen)

regarding g/day, percentage of the substance, percentage of energy...



As outlined in chapter "Indirect Calorimetry", there are various calculation formulas the user can generate via the software settings.

# Basics and Diagnostics

## Threshold Determination

An important aspect of cardiopulmonary exercise testing is the determination of the different thresholds in order to identify, for example the anaerobic threshold, the respiratory compensation point or the Steady State. The respective thresholds are usually displayed as break points in the CPET graphical displays. Unfortunately, ambiguous terms and abbreviations exist in literature which can lead to confusion, misunderstanding (Binder (2008)) or even an incorrect interpretation.

As an example, Wasserman (2012) refers to the first break point as the anaerobic threshold (AT), whereas it is called the aerobic threshold (AE or AeT) in the field of sports medicine / science (for example Kindermann (2004)).

The same applies to the second break point: Wasserman describes it as the Respiratory Compensation Point (RCP), sports medicine, however, usually refers to it as the anaerobic threshold (AT, IAT or AnT).

The classification according to Weber (1997) was also referred to as the anaerobic threshold but was defined to be the point at which RER = 1. This point is located near the **Maximal Lactate Steady State** (MLSS) (or RCP, according to Wasserman) as well rather than near the anaerobic threshold according to Wasserman (AT).

In sports medicine, lactate is often used to identify the different aerobic and anaerobic areas, since it is easy to determine. This offers the ability to quickly evaluate training success (even in the field) and to compare laboratory measurements with field tests by combining CPET and lactate determination.

Figure 1 schematically shows the classification of the lactate and ventilatory thresholds.

Unfortunately, confusion about the terminology is present here as well: Wasserman defines the beginning of the lactate increase as the so-called lactate threshold LT (and as anaerobic threshold), whereas e.g. Mader (1976) refers to the anaerobic threshold with a lactate value of greater than 4 mmol. This has been improved upon with methods such as the determination of MLSS or further anaerobic lactate thresholds (e.g. Heck (1985), Dickhut (1991), Stegmann (1981), Pokan (2004) and others). Accordingly, LT is located near the first break point (VT1 or AT, according to Wasserman), MLSS near the second break point (VT2 or RCP, according to Wasserman).

Unfortunately, neither the second nor the third (if reached) break point can be recognized directly in the lactate curve.

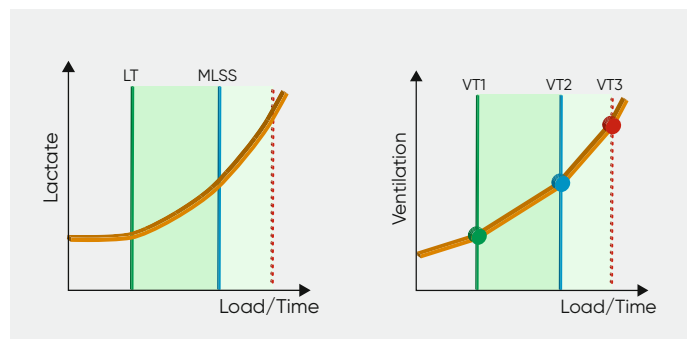


Fig. 1: Connection between lactate and ventilatory thresholds

Training ranges are usually specified by means of VT2 (AT according to sports medicine, RCP according to Wasserman). However because of the confusion mentioned above some therapists / trainers accidentally use the anaerobic threshold according to Wasserman for the creation of a training schedule which are consequently at too low intensity and effectively useless.

In 2012, a major CPET working group (Westhoff 2013) decided to refer to the thresholds as VT1 and VT2 according to their emergent order rather than continuing to name them aerobic, anaerobic or respiratory compensation point in order to avoid this confusion in the future. These terms seem to have gained wider international acceptance which prompted us to implement these definitions into all our current software versions.

Additionally, it is possible to recognize a third break point (VT3), primarily in high-performance athletes. So far, this point is hardly described in literature and its meaning is still not entirely clear (see some notes further below). Nevertheless we implemented this point as VT3 in our software. Sport medicine / science occasionally refers to this point as "Respiratory Compensation Point", "Panic breathing" or as "Hot ventilation".



Following table shows a summary of the various terms and how they are implemented in our programs:

SeS, JLAB	JLAB <V 5.72	Binder	Wasserman	Sports medicine	Lactate
VT1	AE	AE	AT	AeS, AeT,	1. Lactate threshold (Beginning of the lactate to increase), LT
VT2	AT	AT	RCP	AnT, IAS, MLSS, IAT, ...	2. Lactate threshold, MLSS (ca. 4 mmol) or various other methods such as Dickhut, Stegmann etc.
VT3	RCP	??	??	Sometimes referred to as Hot Ventilation (HV), sometimes labelled RCP	–

**Fig. 2: Summary of the various threshold terms**

The respective abbreviations stand for:

VT1, VT2, VT3	=	First, second and third ventilatory threshold
AE, AeS, AeT	=	Aerobic threshold
AT, AnT	=	Anaerobic threshold
IAS, IAT	=	Individual anaerobic threshold
LT	=	Lactate threshold (beginning of the lactate to increase)
MLSS	=	Maximum Lactate Steady State
RCP	=	Respiratory Compensation Point
HV	=	Hot Ventilation ("extreme" hyperventilation)

As already mentioned and as it can be seen in the table above, different and partially misleading terms are used for the lactate thresholds as well. It is therefore recommended to use the term (analogically for the ventilatory thresholds) LT1 for the beginning and LT2 for the second threshold. Furthermore, it should be noted which of the more than 60 different concepts have been used for the lactate threshold determination.

With Cardiopulmonary Exercise Testing it is possible to recognize the respective thresholds by means of more or less distinctive break points. Depending on the selected parameters they can vary from: obvious; to ambiguous; to not apparent, in the various graphical displays. For this reason, it is recommended to observe several plots simultaneously. For clarity, it should be noted that the "thresholds" we usually refer to are in fact transition areas which do not appear as an individual data point but appear as a transition area due to the body's different control mechanisms. Therefore, it is possible that these points can seem to occur in slightly different positions from one graph to another. The averaging of the individual parameter is significant as well (for example: over 8 breaths; 30 seconds or no averaging at all). According to our experience, a moving average over 8 to 10 breaths has proven its worth.

Details of the appearance of the respective break points shall be given in the following.

# Basics and Diagnostics

## Ventilatory Threshold VT1

Typically, we burn a mixture of approx. 50% fat and 50% carbohydrates during rest and low intensity exercise. (As proteins only contribute a little to energy production, they will be disregarded for the purpose of this discussion). Metabolic oxidation of fats produces approx. 700 mL CO<sub>2</sub> from 1000 mL oxygen. Thus, the RER (=V'CO<sub>2</sub>/V'O<sub>2</sub>) is ~0.70. The breakdown of carbohydrates results in 1000 mL of CO<sub>2</sub> being produced of 1000 mL O<sub>2</sub> resulting in an RER of ~1.00. The average RER with a mixed metabolism is consequently approx. (0.70 + 1.0) / 2 = 0.85. The V-Slope graphic V'CO<sub>2</sub>(V'O<sub>2</sub>) – Fig. 3 therefore shows a linear increase with a slope of less than 1 in the lower area (approx. 0.85). If the same amount of CO<sub>2</sub> is produced as oxygen was utilized the RER would be 1.00, as represented by the dashed line.

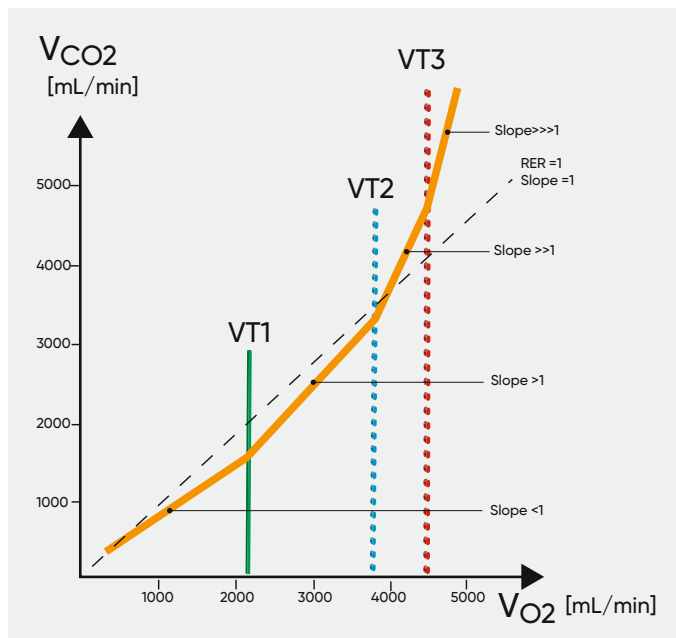


Fig. 3: V-Slope Graphic illustrating the emergence of the break points

An increased load results in an increased energy requirement. The body recognizes the demand for a more efficient performance so it tries to metabolize more carbohydrates and less fat. As a consequence, more CO<sub>2</sub> per oxygen uptake is produced. At the same time, the anaerobic glycolysis is initiated which causes CO<sub>2</sub> (and lactate) to be released as well. As a result, the curve is slowly approaching the RER=1 line and the slope is > 1 from that moment on. The resulting increase in CO<sub>2</sub> production drives a proportionate increase in ventilation due to the causal relationship between these two parameters.

Thus, the break point VT1 is also evident in several other graphics (Fig. 4), such as

- V'CO<sub>2</sub> (V'O<sub>2</sub>) (V-Slope graphic)
- EqO<sub>2</sub> (Time, Load)
- V'E (Time, Load, V'O<sub>2</sub>)
- V'CO<sub>2</sub> (Time, Load, V'O<sub>2</sub>)
- PETO<sub>2</sub> (Time, Load)

For verification, it is often helpful to consult further parameters.

- EqCO<sub>2</sub> (Time, Load) – decreases and subsequently turns into a plateau
- PETCO<sub>2</sub> (Time, Load) – increases and subsequently turns into a plateau

As previously mentioned, anaerobic metabolism starts to contribute a larger proportion of energy causing increased production of CO<sub>2</sub> and lactate. However, the CO<sub>2</sub> and lactate amount is still so small that the body is able to metabolize it (if the load does not increase any further).

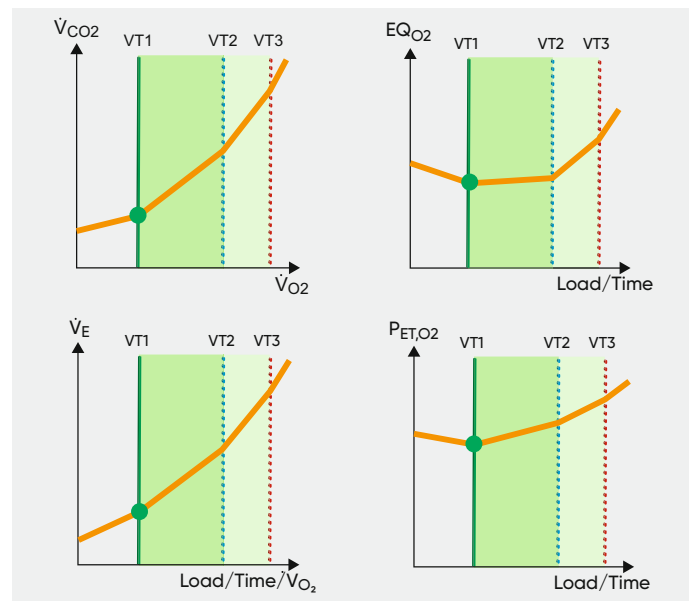


Fig. 4: Graphics for the determination of VT1

## Ventilatory Threshold VT2

When the subject is no longer able to provide the muscles with enough oxygen for the generation of energy and the (still relatively low) additional anaerobic metabolism does not provide enough energy any more, the body will intensify the anaerobic metabolism. Thus, further lactate and, with it, further CO<sub>2</sub> is produced. This transition becomes visible in form of another break point (VT2, Fig. 3 and Fig. 5), since the body will disproportionately increase ventilation due to the increasing metabolic acidosis:

- $\dot{V}'E$  ( $\dot{V}'CO_2$ )
- $EqCO_2$  (Time, Load)
- $\dot{V}'E$  (Time, Load,  $\dot{V}'O_2$ )
- $\dot{V}'CO_2$  (Time, Load,  $\dot{V}'O_2$ )
- $PETCO_2$  (Time, Load)

Further parameters can be used for verification:

- $EqO_2$  (Time, Load) – further break point upwards
- $PETO_2$  (Time, Load) – further break point upwards

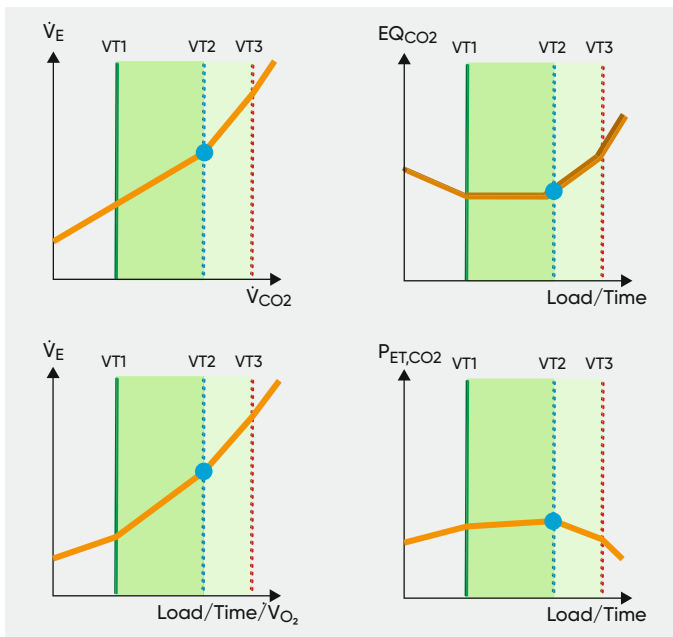


Fig. 5: Graphics for the determination of VT2

## Ventilatory Threshold VT3

If the subject is able to continue to exercise far beyond VT2, a third break point (VT3, Fig. 3) can be observed. This threshold, however, can only rarely ever be reached: Due to the continually increasing acidosis, the subject is no longer able to control his breathing and starts to hyperventilate (sometimes also referred to "panic breathing" or "hot ventilation"). Now, the breathing only aims to compensate the metabolic acidosis and to eliminate the accumulated  $CO_2$  as quickly as possible. This can often be noticed by the breathing sound and can be recognized in different graphics (Fig. 6):

- $\dot{V}'E$  ( $\dot{V}'CO_2$ )
- $EqO_2$ ,  $EqCO_2$  (Time, Load)
- $PETO_2$ ,  $PETCO_2$  (Time, Load)
- $\dot{V}'E$  (Time, Load,  $\dot{V}'O_2$ ) or  $\dot{V}'CO_2$  (Time, Load,  $\dot{V}'O_2$ )

However, the usage of this third break point regarding further evaluation, interpretation or exercise prescription is still investigated and will be addressed in the future.

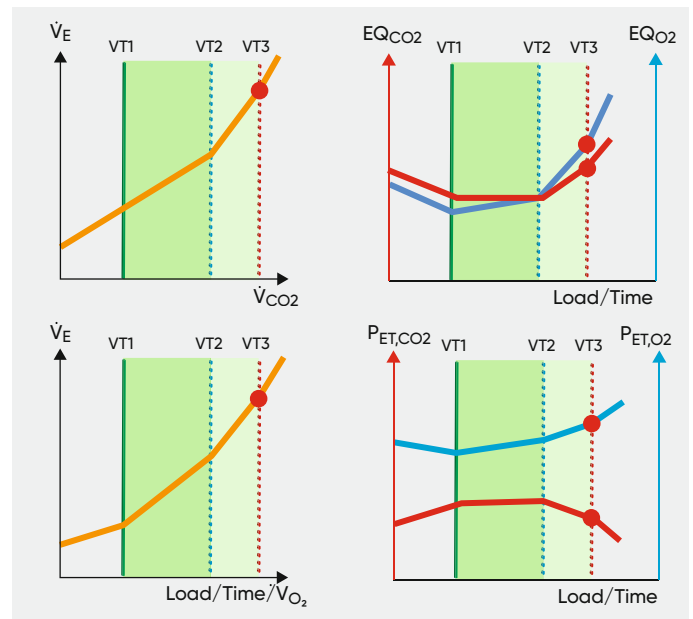


Fig. 6: Graphics for the determination of VT3

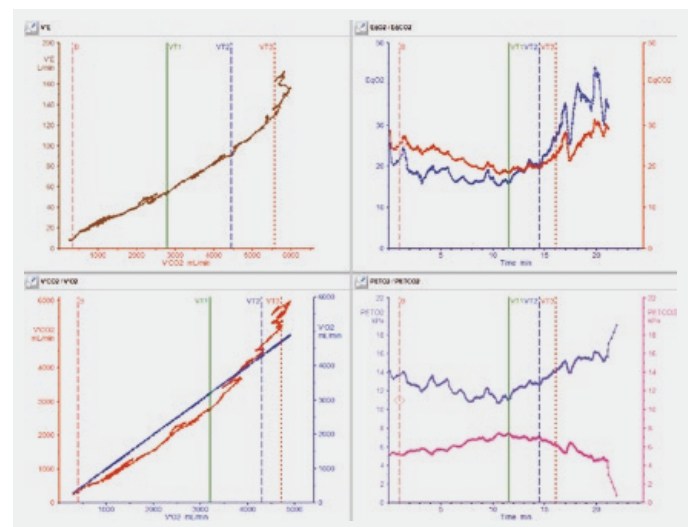


Fig. 7: Measurement example for all three thresholds

## Comparison between Ventilatory and Lactate Thresholds

The thresholds VT1 and VT2 mentioned above can be assigned to the respective lactate thresholds LT and MLSS provided that the measurements were performed under identical conditions, e.g. identical load profile, identical load device. With such comparisons, however, it needs to be observed that

- measurements including lactate determination are usually performed with a step protocol
- CPET measurements are usually performed with a ramp protocol, causing  $\dot{V}'E$ ,  $\dot{V}'O_2$ ,  $\dot{V}'CO_2$  and further parameters to follow the load

# Basics and Diagnostics

- the lactate determination is often interrupted in order to take blood samples. During this interruption period, a short “rest” occurs ( $\dot{V}E$ ,  $\dot{V}O_2$ , HR and other parameters decrease) and the lactate will partially be metabolized.
- when resuming the exercise, the subject must provide a large amount of anaerobic energy at first in order to reach the steady state (or quasi-steady-state) again.

## Threshold Application

As already outlined, the respective threshold can be used to determine various ranges such as aerobic exercise (below VT1) or the different aerobic-anaerobic transitions (VT1 and VT2). A reduced VT1 with unhealthy subjects can be used for further differential diagnosis. In particular, the area between VT1 and VT2 or VT2 itself can serve exercise prescription purposes (e.g. Kindermann (2004)).

As the length of this chapter is limited, further information on this topic is available in the respective literature (e.g. Jones (1999), Cooper (2001), Pokan (2004), Rühle (2008) Wasserman (2012), Kroidl (2013)).

Furthermore, conclusions regarding the potential endurance capacity can be made from the different VT's: This capacity is important in the context of occupational medicine, retiring or end of working life-span. Comparing these values with tables on task specific exercises can inform decisions for retirement or occupational redeployment. However it should be noted, that - similar to the thresholds - misunderstandings regarding the term “endurance capacity” exist:

Some publications (especially in sports medicine) refer to MLSS (consequently, VT2 in CPET) as the endurance limit whereas other publications choose VT1 as the limitation (e.g. Hollmann (1959) – point of best efficiency).

Although it is obvious that exercising at the MLSS without completely exploiting the anaerobic reserves is possible for a specific amount of time (approx. 30–60 minutes), a subject cannot be expected to be able to continue this throughout an entire day.

Unfortunately, various policies (for example retirement) do not mention the thresholds at all but focus on  $\dot{V}O_2\max$ , or only discuss the “anaerobic threshold”. The question of which threshold (VT1 or VT2) is actually meant often goes unanswered.

In my opinion, the steady state limit (for an 8-hours-day) is most closely approximated by VT1, whereas it is located near VT2 with regards to short term exercise (less than 1 hour).

### References:

- **Binder R.K., Wonisch M., Corra U., Cohen-Solal A., Vanhees L., Saner H., Schmid J.-P.:** Methodological approach to the first and second lactate threshold in incremental cardiopulmonary exercise testing. *Eur J of Cardiovascular Prevention and Rehabilitation* 6\_15 (2008); 726–734.
- **Cooper C.B., Storer T.W.:** Exercise testing and interpretation. Cambridge University Press (2001). ISBN: 0-521-64842-4.
- **Dickhuth H.-H., Huonker M., Munzel T. et al.:** Individual anaerobic threshold for evaluation of competitive athletes and patients with left ventricular dysfunctions. In: Bachl N, Graham TE, Löllgen H (Hrsg.): *Advances in ergometry*. Berlin, Heidelberg, New York, Springer (1991); 173–179.
- **Heck H., Mader A., Hess G., Mücke S., Müller R., Hollmann W.:** Justification of the 4–mmol/l lactate threshold. *Int J Sports Med* 6 (1985); 117–130.
- **Hollmann W.:** The relationship between pH, lactic acid, potassium in the arterial and venous blood, the ventilation, PoW and puls frequency during increasing spirometric work in endurance trained and untrained persons. Chicago, 3rd Pan-American Congress for Sports Medicine (29.11.1959).
- **Jones N.L.:** *Clinical Exercise Testing*. 4th Edition, W.B. Saunders Company (1997). ISBN: 0-7216-6511-x.
- **Kindermann W.:** Anaerobe Schwelle. Salzburg, Austria (1986) Bd.-Hrsg. ISBN: 3886032388. *Dt Zeitschrift für Sportmedizin* 55 (2004); 161–162.
- **Kroidl R.F., Schwarz S., Lehnigk B., Fritsch J.:** *Kursbuch Spiroergometrie – Technik und Befundung verständlich gemacht*. 3. Auflage Thieme Verlag (2014). ISBN: 978-3-13-143443-2.
- **Mader A, Liesen H, Heck H., Philippi H., Rost R., Schuerch P., Hollmann W.:** Zur Beurteilung der sportartspezifischen Ausdauerleistungsfähigkeit im Labor. *Dtsch Z Sportmed* 27 (1976); 80–112.
- **Pokan R., Förster H., Hofmann, P., Hörtnagl H., Ledl-Kurkowski, Wonisch M.:** *Kompendium der Sportmedizin*. SpringerWienNewYork (2004). ISBN: 3-211-21235-1
- **Rühle K.-H.:** *Praxisleitfaden der Spiroergometrie*. 2. überarbeitete und erweiterte Auflage. Kohlhammerverlag (2008). ISBN: 978-3-17-018053-6.
- **Stegmann H., Kindermann W., Schnabel A.:** Lactate kinetics and individual anaerobic threshold. *Int J Sports Med* 2 (1981) 160–165.
- **Wasserman K., Hansen J.E., Sue D.Y., Stringer W.W., Sietsema K.E., Sun X-G., Whipp B.J.:** *Principles of Exercise Testing and Interpretation*. 5th edition (2012). Lippincott Williams & Wilkins. ISBN-13: 978-1-60913-899-8.
- **Weber K.T.:** What can we learn from exercise testing beyond the detection of myocardial ischemia. *Clin Cardiol* 20 (1997); 684–696.
- **Westhoff M., Rühle K.-H., Greiwing A., Schomaker R., Eschenbacher H., Siepmann M., Lehnigk B.:** Ventilatorische und metabolische (Laktat-)Schwellen. *Dtsch Med Wochenschrift* 138 (2013); 275–80.


## Indirect Calorimetry

### Basics

Parameters measured or calculated with cardiopulmonary exercise testing ( $\dot{V}O_2$ ,  $\dot{V}CO_2$  as well as RER) and the additional determination of urinary nitrogen (if applicable), can be used to obtain information on the energy expenditure and the relative contribution of the individual substrates (carbohydrates, fat and proteins). Calculations are based on the respective chemical equations for the oxidation of the various substrates. However, it must be noted that there are different substances (e.g. glucose, disaccharides or polysaccharides for carbohydrates). Consequently, oxygen consumed and carbon dioxide produced varies according to the basic substance. As a result, the energy production as well as the RER can be slightly different depending on the substrates.

Different examples of substances and chemical equations are listed below:


#### 1. Carbohydrates:

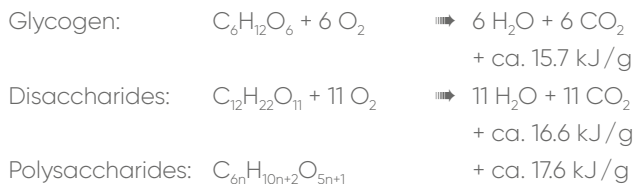

In this case, a specific amount of energy is released per mole with the amount of  $O_2$  required and  $CO_2$  produced being identical (RER = 1.00):


$$RER = 6/6 = 11/11 = \dots = 1.00$$

#### 2. Fats:

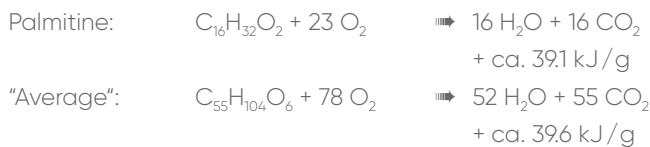

This results in an average RER of about 0.70  
( $16/23 = 0.696$ ;  $55/78 = 0.705$ ).

#### 3. Proteins:

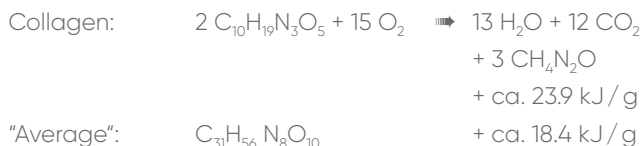


This results in an RER of about 0.80 (=12/15).

The protein metabolism – even though it contributes only little to energy production and is consequently often ignored – can be derived from the urinal nitrogen value:

Approximately 16% of all converted proteins is excreted via the urinal nitrogen which allows for the calculation of the total protein oxidation rate (dP):

$$dP = 6.25 \cdot UN$$

Furthermore, the chemical equations listed above can be used to calculate the amount of oxygen required as well as carbon dioxide produced per mole (or g), e.g. for

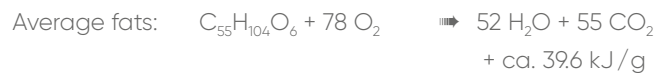

1 mole  $C_{55}H_{104}O_6$  weighs about 860 g.

78 mole  $O_2$  weigh about  $78 \cdot 32 \text{ g}$  and occupy about  $78 \cdot 22.4 \text{ L}$ .

55 mole  $CO_2$  weigh about  $55 \cdot 44 \text{ g}$  and occupy about  $55 \cdot 22.4 \text{ L}$ .

Thus, in order to metabolize 860 g fat, about 1747 L  $O_2$  and 1232 L  $CO_2$  are required. To break down 1 g fat, respectively, about 2.03 L  $O_2$  and 1.43 L  $CO_2$  are required.

As already mentioned, carbohydrates, fats as well as proteins consist of different molecules. Consequently, the metabolization is the respective average of the single molecules.

Depending on the mixture, slightly different values for the caloric value as well as the respective gas proportion ( $\dot{V}O_2$  and  $\dot{V}CO_2$ ) therefore exist.

Substrate	Gas volume required for 1 g of the substrate			Caloric Value	
	$O_2$ [L]	$CO_2$ [L]	RER	[kJ/g]	[kcal/g]
CHO	0.829	0.829	1.000	17.50	4.18
FAT	2.019	1.427	0.707	39.61	9.46
Protein	0.966	0.782	0.810	18.09	4.32
UN	6.040	4.890	0.810	113.05	27.00

#### Average indirect calorimetry values according to Takala (1989)

According to those average values of Takala (1989), the total amount of  $\dot{V}O_2$  or  $\dot{V}CO_2$ , respectively, is calculated as follows:

$$\begin{aligned} \dot{V}O_2 &= 0.829 \cdot CHO + 2.019 \cdot FAT + 0.966 \cdot Prot \\ &= 0.829 \cdot CHO + 2.019 \cdot FAT + 6.040 \cdot UN \end{aligned} \quad (1)$$

$$\begin{aligned} \dot{V}CO_2 &= 0.829 \cdot CHO + 1.427 \cdot FAT + 0.782 \cdot Prot \\ &= 0.829 \cdot CHO + 1.427 \cdot FAT + 4.890 \cdot UN \end{aligned} \quad (2)$$

# Basics and Diagnostics

Solving both equations for CHO, FAT and Proteins provides the following:

$$\text{CHO [g/day]} = 4.12 * V'\text{CO}_2 - 2.91 * V'\text{O}_2 - 2.54 * \text{UN} \quad (3)$$

$$\text{FAT [g/day]} = 1.69 * V'\text{O}_2 - 1.69 * V'\text{CO}_2 - 1.94 * \text{UN} \quad (4)$$

$$\text{Prot [g/day]} = 6.25 * \text{UN} \quad (5)$$

The total energy results of

$$\text{EE[kcal/d]} = 4.18 * \text{CHO} + 9.46 * \text{FAT} + 27 * \text{UN} \quad (6)$$

Using the individual contingents (3, 4, 5), (6) becomes:

$$\begin{aligned} \text{EE[kcal/d]} &= 4.18 * (4.12 * V'\text{CO}_2 - 2.91 * V'\text{O}_2 - 2.54 * \text{UN}) \\ &+ 9.46 * (1.69 * V'\text{O}_2 - 1.69 * V'\text{CO}_2 - 1.94 * \text{UN}) \\ &+ 27 * \text{UN} \end{aligned}$$

$$\text{EE[kcal/d]} = 3.82 * V'\text{O}_2 + 1.22 * V'\text{CO}_2 - 1.99 * \text{UN} \quad (7)$$

Correspondingly, the single energies for CHO, FAT and Proteins can be calculated.

As already mentioned, literature partially uses different caloric values resulting in slight differences in calculations. Other formulas frequently used are those from Frayn (1983) or e.g. those from de V. Weir (1949):

$$\text{EE [kcal/d]} = 3.94 * V'\text{O}_2 + 1.11 * V'\text{CO}_2 - 2.17 * \text{UN}$$

## Note 1:

With the formulas derived above,  $V'\text{O}_2$  and  $V'\text{CO}_2$  are used in L/day.

For the conversion into the units used in cardiopulmonary exercise testing (mL/min), the factor results in  $24 * 60 / 1000 = 1.44$ . As an example, the de Weir formula becomes:

$$\begin{aligned} \text{EE [kcal/d]} &= 3.94 * V'\text{O}_2 + 1.11 * V'\text{CO}_2 - 2.17 * \text{UN} \\ & \quad (V'\text{O}_2, V'\text{CO}_2 \text{ in L/day; UN in g/day}) \end{aligned}$$

$$\begin{aligned} \text{EE [kcal/d]} &= 5.67 * V'\text{O}_2 + 1.60 * V'\text{CO}_2 - 2.17 * \text{UN} \\ & \quad (V'\text{O}_2, V'\text{CO}_2 \text{ in mL/min; UN in g/day}) \end{aligned}$$

## Note 2:

Those calculations for CHO, FAT and proteins are actually intended for steady state measurements only. Consequently, they cannot be used for  $\text{RER} < 0.70$  and  $\text{RER} > 1!$

## Fields of Application

Indirect calorimetry is used in various fields of application. Therefore, it can only be briefly discussed in the context of this special edition.

For more detailed information please refer to the appropriate literature (e.g. Ferrannini (1988), De Lorenzo (2001), AARC (2004), Schols (2014)).


The proper relation between food intake and energy expenditure is essential for a proper, balanced nutrition. By means of indirect calorimetry, the energy expenditure can be determined to adjust the nutrition according to the individual goal (weight gain, weight loss, weight maintenance, or activating the fat burning process, for example).



## 1. BxB (or mixing chamber):

As all required units are available in BxB-Mode as well as with the mixing chamber, a standard BxB measurement already provides the respective parameters.

Thus, the energy production can already be graphically displayed during an exercise test and, as an example, the optimum load for the fat burning process can be determined (Note: the subject needs at least to be in a "Quasi Steady State". A ramp profile is the best approach for the patient to achieve this state.)

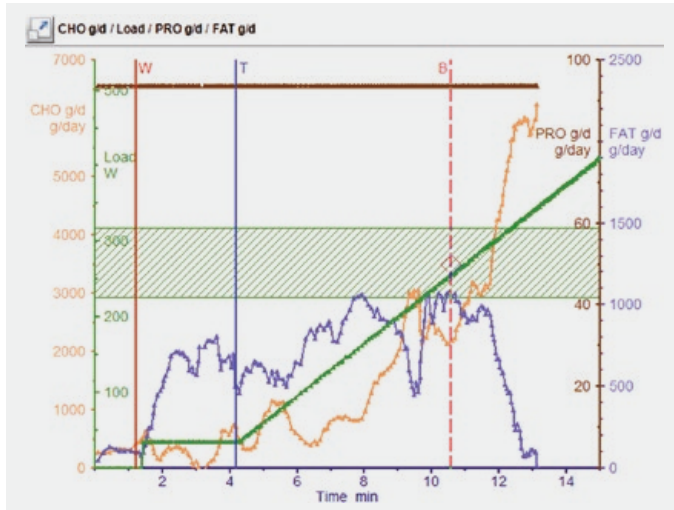

Indirect Calorimetry during exercise

## 2. Canopy (also see "Option Canopy"):

The canopy dilution method is a patient friendly test procedure and the gold standard to determine the resting energy expenditure: The subject lies for ca. 20–30 min while relaxing (or even sleeping) under a canopy. The gas exhaled is analyzed for  $\dot{V}O_2$ ,  $\dot{V}CO_2$ , RER... via a blower system and is evaluated regarding resting energy expenditure, carbohydrates, fats and proteins according to the formulas mentioned above.



Indirect calorimetry using the canopy dilution method

## 3. Patients with additional oxygen supply:

Indirect calorimetry is of great benefit especially for ventilated patients:

As a patient's healing process is delayed in case of both supernutrition as well as malnutrition, a detailed knowledge of the resting energy expenditure is necessary, especially as the ventilator carries out a part of the work of breathing, which is unknown.

By means of indirect calorimetry, the ventilated patient's resting energy expenditure can be determined as well. Furthermore, energy production can be divided according to fats, carbohydrates and proteins in order to optimally adjust enteral nutrition.

Due to the large variety of ventilators and ventilation modes (for example bias flow with increased oxygen concentration), the adaptation to a ventilator may differ from case to case. Especially with an increased oxygen concentration, the Haldane transformation turns out to be questionable (Eschenbacher (2016)).

### References:

- **AARC:** Clinical Practice Guideline - Metabolic Measurement Using Indirect Calorimetry During Mechanical Ventilation - Revision and Update. Respiratory Care 9\_49 (2004).pdf
- **De Lorenzo A., Tagliabue A., Andreoli A., Testolin G., Comelli M., Deurenberg P.:** Measured and predicted resting metabolic rate in Italian males and females, aged 18 ± 59 y. Eur Journal of Clin Nutr 55 (2001); 208-214.
- **de V. Weir J.B.:** New methods for calculating metabolic rate with special reference to protein metabolism. J Physiol. 109 (1949); 1-9.
- **Eschenbacher H.:** Haldane and Eschenbacher transformation. White Paper RD8168 (0118/PDF). Vyaire (2018)
- **Ferrannini E.:** The Theoretical Bases of Indirect Calorimetry: A Review. Metabolism 3\_37 (1988); 287-301.
- **Frayn K.N.:** Calculation of substrate oxidation rates in vivo from gaseous exchange. J Appl Physiol 2\_55 (1983); 628-634.
- **Schols A.M., Ferreira I.M., Franssen F.M., Gosker H.R., Janssens W., Muscaritoli M., Pison C., Rutten-van Mölken M., Slinde F., Steiner M.C., Tkacova R., Singh S.J.:** Predicting basal metabolic, new standards and review of previous work. Eur Respir J 44 (2014); 1504-1520.
- **Takala J., Meriläinen P.:** Handbook of Gas Exchange and Indirect Calorimetry. Libris Oy, Datex Finland (1989); Doc No. 876710.

## Haldane and Eschenbacher Transformation\*

### Background and History

Over 100 years ago, Haldane (1912) proposed the Haldane transformation (HT) as the calculation basis for cardiopulmonary exercise tests which has been used ever since. However, notes can repeatedly be found in literature at least questioning the validity of this calculation. Prieur (2002), for example, reports that increased oxygen uptakes with hyperoxia measurements can only partially be explained. Whereas Stanek (1979) describes that the measurement of oxygen uptake with elevated oxygen concentrations – calculated by means of the HT – is incomprehensively increased whereas calculations via blood gas analysis and the direct Fick method, however, do not confirm this. In 1986 I attended a congress, "Methodische Fragen zur Indirekten Kalorimetrie", in Austria, where the methods of indirect calorimetry measurements were discussed (Kleinberger (1986)).

Of particular interest: Why during cardiopulmonary exercise testing measurements seem accurate with normal breathing but are implausible with elevated  $FIO_2$  concentrations?

The group reached following conclusions:

- When  $FIO_2 < 40\%$ , measurements seem accurate.
- With  $FIO_2$  between 40% and 60%, a careful calibration is required to achieve results that are plausible, although not always.
- When  $FIO_2$  is between 60% and 80%, most values are not plausible.
- For  $FIO_2$  above 80%, all values are implausible.
- At  $FIO_2$  of 100%, no calculation of  $VO_2$  is possible.

Similar conclusions can be found in the "Handbook of Gas Exchange and Indirect Calorimetry" published by the Finnish company Datex (Takala (1989)).

During the same period, in the late 1980s, one manufacturer of cardiopulmonary exercise testing even withdrew the system from the market due to similar concerns regarding results that were either implausible and/or not reproducible.

In 1987 JAEGER™ received similar complaints from customers in Italy and South Africa that the values delivered from the EOS-Sprint were implausible at elevated  $FIO_2$  concentrations.

I carefully repeated the tests and got the same results as those reported from Italy and South Africa. I discovered that these inaccurate results seemed to be a general problem with the Haldane transformation. I set about to solve this problem by creating a new set of formulas.

These new formulas delivered plausible results over the whole range, even at  $FIO_2 = 100\%$ .

**Note:** The inspired and expired volumes ( $V'I$ ,  $V'E$ ) are expressed in BTPS, while  $V'O_2$  and  $V'CO_2$  are expressed in STPB. For simplification of the formulas, the conversion factors are ignored in the following discussions. The change in the water vapor content between inspiration and expiration can be ignored as the analyzed gases are conditioned (dried) before analysis. Furthermore,  $FICO_2$  (normally ca. 0.03 – 0.05%) is ignored as well.

### The Haldane Transformation (HT)

(see e.g. Consolazi (1963))

Oxygen uptake ( $V'O_2$ ), carbon dioxide output ( $V'CO_2$ ) as well as nitrogen exchange ( $V'N_2$ ) are calculated as the difference between inspired and expired volumes. The following basic calculations are used:

$$V'O_2 = FIO_2 \times V'I - FEO_2 \times V'E \quad (1)$$

$$V'CO_2 = FECO_2 \times V'E - FICO_2 \times V'I \quad (2)$$

$$V'N_2 = FIN_2 \times V'I - FEN_2 \times V'E \quad (3)$$

with:

FI = mean inspired gas fractions of  $O_2$ ,  $CO_2$  and  $N_2$

FE = mean expired gas fractions of  $O_2$ ,  $CO_2$  and  $N_2$

$V'I$  = inspired volume

$V'E$  = expired volume

\* The term "Eschenbacher Transformation" was coined by our Italian representative after verifying that the formula I developed was giving plausible values over the whole range of  $FIO_2$ , even at 100%.



During ergospirometry, traditionally only the expired volume is measured, while the inspired volume is calculated via the Haldane transformation. Haldane made the assumption that there is no nitrogen exchange:

$$V'N_2 = 0 \quad (4)$$

With this assumption, equation (3) leads to:

$$V'I = V'E \times (FEN_2 / FIN_2) \quad (5)$$

The low concentration gases in the air (e.g., helium or argon) act like nitrogen and can be neglected or added to the nitrogen content.

This results in the following two equations:

$$FIN_2 + FIO_2 + FICO_2 = 1 \quad (6)$$

$$FEN_2 + FEO_2 + FECO_2 = 1 \quad (7)$$

or:

$$FIN_2 = 1 - FIO_2 - FICO_2 \quad (6a)$$

$$FEN_2 = 1 - FEO_2 - FECO_2 \quad (7a)$$

(6a) and (7a) in (5) leads to the following result:

$$V'I = V'E \times (1 - FEO_2 - FECO_2) / (1 - FIO_2 - FICO_2) \quad (8)$$

(8) in (1) therefore gives:

$$V'O_2 = V'E \times kH \times FIO_2 - V'E \times FEO_2 \quad (9)$$

with the Haldane correction factor:

$$kH = (1 - FEO_2 - FECO_2) / (1 - FIO_2 - FICO_2) \quad (10)$$

## Discussion of the Haldane transformation

Under normal conditions we can expect for a constant workload, that below the ventilatory threshold 2 (VT2) the same oxygen uptake is needed as well as the same carbon dioxide is produced, independent of the inspired  $FIO_2$ .

So the difference of the gas fractions should be constant:

$$DFO_2 = FIO_2 - FEO_2$$

respectively

$$DFCO_2 = FECO_2 \text{ (with } FICO_2 = 0)$$


Example: At a load of 40 W, the following measuring results are expected:

$$V'E = 20 \text{ L/min}$$

$$DFCO_2 = 4\%$$

$$DFO_2 = 4.8\%$$

Using the Haldane transformation with these values and a varying  $FIO_2$ , we will get the following results:

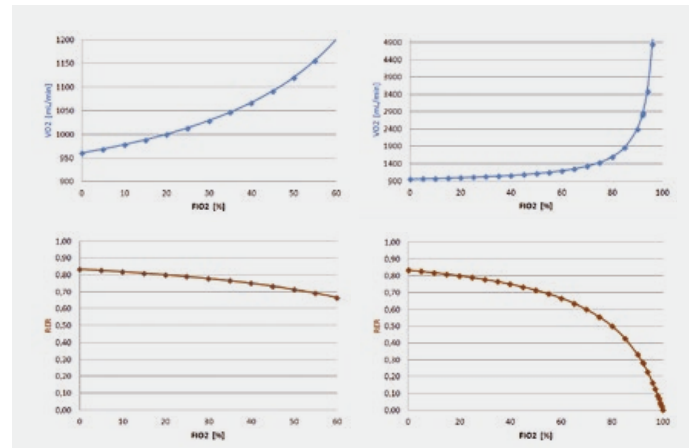

$V'O_2$  and RER as a function of  $FIO_2$  with the Haldane transformation for a typical 40 W exercise. Note the scaling. **Left side:**  $V'O_2$  increases from ca. 1000 mL/min at 20%  $FIO_2$  to around 1200 mL/min at 60%  $FIO_2$ , while RER decreases from 0.80 down to <0.70. **Right side:**  $V'O_2$  increases to > 5000 mL/min at 97%  $FIO_2$  and goes to infinity, while RER goes down to zero.

The following conclusions can be derived:

- With  $FIO_2$  approaching 100% the calculated  $V'O_2$  goes to infinity. This is due to the Haldane transformation, which is obviously not valid. That also seems to be the reason why in "Principles of Exercise Testing and Interpretation" (Wasserman (2012)) all cases with oxygen breathing do not show any data for  $V'O_2$ , RER and other depending parameters.
- For  $FIO_2$  going to 0% the oxygen uptake gives the same value as if  $V'I = V'E$ . This, however, is for example with  $V'O_2 = 960 \text{ mL/min}$   
 $V'CO_2 = 800 \text{ mL/min}$   
in contradiction to:  
 $V'I - V'E = V'O_2 - V'CO_2 = 160 \text{ mL/min}$  or  $V'I \neq V'E$
- Due to the Haldane transformation, the  $V'I$  (and therefore also the inspiratory tidal volume  $VTin$ ) should increase dramatically with a high  $FIO_2$ , for example above to  $VTin > 2 \times VTex$  at 99.2%  $FIO_2$ . However, such differences could not be measured and would cause an enormous drift in the spirogram, which could not be observed as well.

# Basics and Diagnostics

- As the formula is neither valid for  $FIO_2$  nearing 0%, nor for  $FIO_2$  going to 100%, and measurements often show, that the results are already questionable at  $FIO_2$  of about 50% (1000 mL/min at 20%, 1125 mL/min and RER = 0.7 at 50%) the question is: For which  $FIO_2$  can the Haldane transformation be applied at all?
- Last but not least, many publications indicate that there is also a nitrogen exchange during the respiration (both a retention as well as a production, is possible depending for example on the content of the last meal and measurement time after the last meal - Wilmore (1973)). This of course is in contradiction to the assumption (4), that  $V'N_2 = 0$ .


## New considerations for Eschenbacher transformation (ET)

Both the implausible values at elevated oxygen breathing as well as the fact that the Haldane transformation cannot be applied at 100% oxygen breathing made it necessary to develop a new calculation which

- is not based on the assumption that  $V'N_2 = 0$ .
- still takes into account that for RER "unequal" to 1,  $V'I$  is different to  $V'E$ .
- calculates plausible values also at elevated  $FIO_2$ .
- even allows to calculate  $V'O_2$  at  $FIO_2 = 100\%$ .


### Measurements at normal room air

At normal ambient conditions ( $FIO_2 = 20.93\%$ ), both calculations deliver the same values within the measurement accuracy.

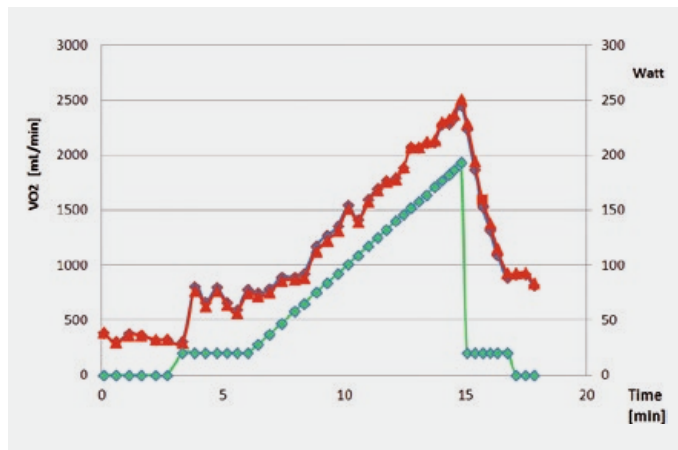

$V'O_2$  at normal conditions ( $FIO_2 = 20.93\%$ ), evaluated with the HT (blue) and ET (red). No significant differences between both calculations can be observed.

Also the Bland-Altman comparison shows a good agreement. In order to facilitate the comparison, the following graphic does not show the difference between both methods but the respective difference to the mean value:

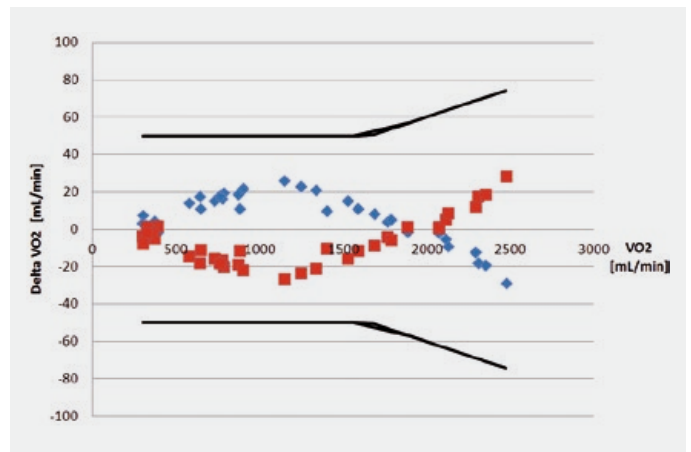


$V'O_2$  deviation from the mean value for HT (blue) and ET (red) at  $FIO_2 = 20.93\%$ . RER = 1 is reached at ca.  $V'O_2 = 2000$  mL/min. While the HT  $V'O_2$  is a bit higher for RER < 1 and lower for RER > 1, the ET is a bit lower for RER < 1 and higher for RER > 1. Both calculations, however, are within the measurement accuracy (solid lines).

The  $V'O_2$  with the HT is a bit higher at RER < 1 and a bit lower at RER > 1, while the ET shows the opposite tendency. Therefore, RER will also show small differences between HT and ET, but both deviations are within the measurement accuracy.

### Measurements at elevated $FIO_2$

More obvious are the differences at higher  $FIO_2$  values for  $V'O_2$  and RER:

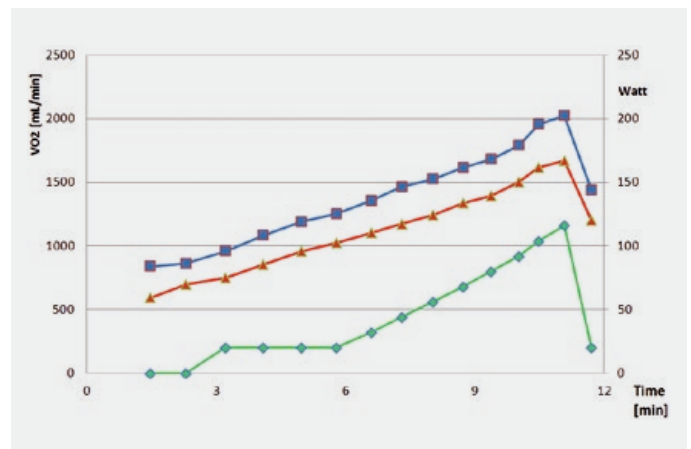


$V'O_2$  measurement at ca.  $FIO_2 = 60\%$ ; while the HT (blue) seems to overestimate  $V'O_2$  (e.g. 1790 mL/min at 90 W) whereas the ET (red) delivers more plausible values (1505 mL/min at 90 W).

According to Wasserman (2012), for  $\dot{V}O_2$  the following is usually expected:

$$\dot{V}O_2 = 5.8 \times BW + 151 + 10.3 \times W$$

which leads in a measurement to a value of ca. 1530 mL/min at 90 W (BW = 79 kg). While the ET calculates a  $\dot{V}O_2$  value close to this expected value (1505 mL/min), the HT seems to overestimate the  $\dot{V}O_2$  at 90 W by ca. 250 mL/min.

The difference between the ET and the HT is even more obvious when comparing the resulting RER:

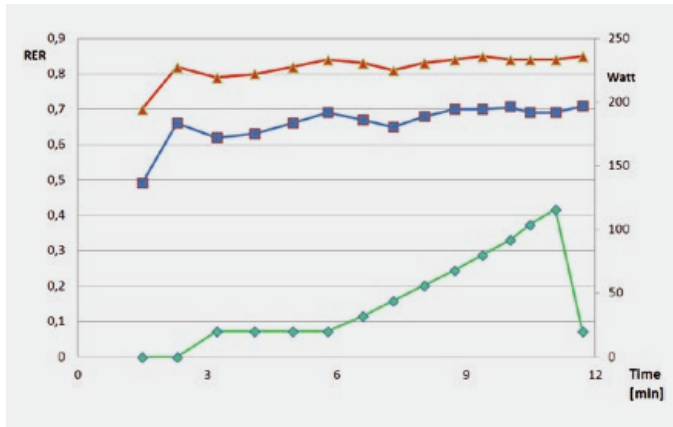

RER measurement at ca.  $FI_{O_2} = 60\%$ ; due to the overestimation of  $\dot{V}O_2$  via the HT (blue), the RER remains implausible below 0.70 almost over the whole measurement, while the ET (red) delivers plausible RER values (between ca. 0.80 and 0.85).

While the ET delivers plausible values for RER (between ca. 0.80 und 0.85), the HT calculates RER values which are unrealistic: The RER with the HT remains below 0.70 (entirely fat burning) even after the wash-in period, which is physiologically impossible.

## Measurements at 100% $FI_{O_2}$ breathing

While the HT does not permit calculating  $\dot{V}O_2$  at 100%  $FI_{O_2}$  at all, the ET delivers plausible values even at 100% oxygen breathing. At the moment, several studies are conducted regarding this issue in order to confirm those aspects and to compare the HT and ET with an elevated  $FI_{O_2}$  regarding their plausibility.

## Conclusion

The HT seems to be limited to  $FI_{O_2}$  values close to room air. Higher  $FI_{O_2}$  values will create significant deviations and HT cannot be used at 100% oxygen breathing.

In contrast, the ET delivers plausible values over the whole  $FI_{O_2}$  range, even when breathing 100% oxygen.

This therefore raises the question, if the Haldane transformation should actually be used due to this limitation, or if it should be skipped from now on?

During my investigations I also had the following impression (though this needs to be investigated in detail, even if it can already be explained by the HT assumption): In the case of nitrogen production, the HT calculation is already implausible at  $FI_{O_2} < 50\%$ , while in case of nitrogen retention the HT seems to deliver more plausible values even at higher  $FI_{O_2}$ .

A change between nitrogen retention and nitrogen production strongly depends on the last meal itself as well as on the time between the meal and the measurement.

Therefore, at least at higher  $FI_{O_2}$  the HT will lead to large fluctuations and non-reproducible results, even with the same patient.

## References:

- **Consolazi CF, Johnson RE, Pecora LJ.:** Physiological Measurements of Metabolic Functions in Man. New York: McGraw-Hill (1963).
- **Haldane JS.:** Methods of air analysis. Charles Griffin & Co Ltd., JB Lippincott Co., Philadelphia (1912).
- **Kleinberger G, Eckart J.:** Methodische Fragen zur indirekten Kalorimetrie. Salzburg, Austria (1986), Bd.-Hrsg. ISBN: 3886032388.
- **Prieur F1, Benoit H, Busso T, Castells J, Geysant A, Denis C.:** Effects of moderate hyperoxia on oxygen consumption during submaximal and maximal exercise. Eur J Appl Physiol. 3\_88 (2002); 235-42.
- **Stanek KA, Nagle FJ, Bisgard GE, Byrnes WC.:** Effect of hyperoxia on oxygen consumption in exercising ponies. J Appl Physiol Respir Environ Exerc Physiol. 6\_46 (1979); 1115-8.
- **Takala J, Meriläinen P.:** Handbook of Gas Exchange and Indirect Calorimetry. Datex Finland (1989); Doc No. 876710.
- **Wasserman K, Hansen JE, Sue DY, et al.:** Principles of Exercise Testing and Interpretation. 5th ed. Lippincott Williams & Wilkins (2012). ISBN-13: 978-1-60913-899-8.
- **Wilmore JH, Costill DL.:** Adequacy of the Haldane transformation in the computation of exercise  $\dot{V}O_2$  in man. J Appl Physiol.(1973);1-35.

# Our CPET History - Examples from 1956-1994



JAEGER  
Ergometer [1956]




SensorMedics  
MMC Horizon [1981]



JAEGER  
Portable [1988]



JAEGER  
Ergo "Glockenspirometer" [1965]

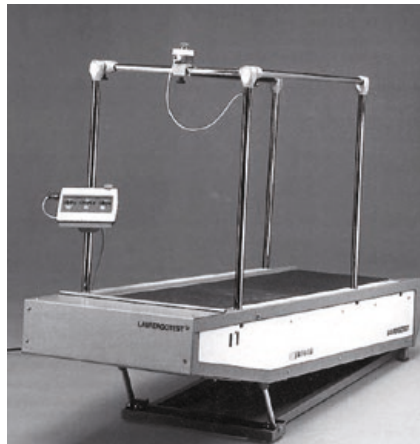

JAEGER  
Laufergotest [1984]



JAEGER  
EOS-Sprint [1989]



Mijnhardt  
Oxycon Portabe [1973]



JAEGER  
Ergo-Oxyscreen [1986]




Mijnhardt  
Oxycon Sigma [1990]



JAEGER  
Ergo-Pneumotest [1976]



SensorMedics  
SMC 2900 [1988]

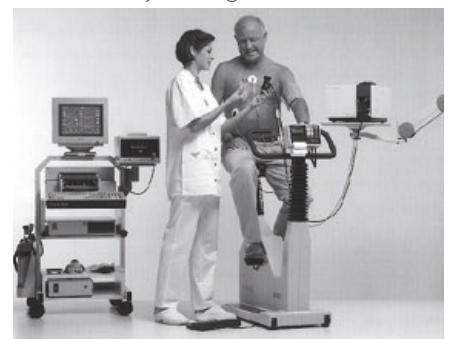

JAEGER  
Oxycon Alpha [1994]

# Examples from 1995–2014



SensorMedics  
Vmax [1995]



JAEGER  
Oxycon Mobile [2002]



JAEGER  
VIASprint [2004]



JAEGER  
LE 200 [1999]



SensorMedics  
Vmax Encore [2004]



CareFusion  
MS-CPX [2008]



JAEGER  
Oxycon Pro [2000]



CareFusion  
Vyntus CPX [2014]



SensorMedics  
Vmax Spectra [2001]

The date of the picture does not necessarily indicate the device's first release.

# Device Presentation

## Vyntus™ WALK



Vyntus WALK is the mobile solution to easily, safely and comfortably perform **6-Minute Walk Test** for assessing the functional exercise capacity of your patients. It can be seamlessly integrated with the diagnostic software platform of SentrySuite for central data management and reporting.

### 6-Minute Walk Test: Small test with a great impact

A 6-minute walk test measures the total distance a patient is able to walk in a period of 6 minutes. This type of exercise testing has consistently been proving value in various medical fields ranging from respiratory, cardiology, physical therapy and rehabilitation medicine to rheumatology, geriatrics and even neurology. The test is especially useful for preoperative and postoperative evaluation, for measuring the response to a therapeutic intervention and for monitoring the medical condition of a patient.\* Additionally, studies show that the 6-minute walk test highly correlates with the maximum oxygen uptake ( $\dot{V}O_2\text{max}$ ) with COPD and cardiac insufficiency and therefore contributes significantly to evaluation. Due to the integration of the latest mobile technology, the 6-minute walk test is quick and easy to perform almost anywhere, making it the preferred test method of both patients and healthcare professionals. High reproducibility, as well as validated and reliable test results and proven diagnostic and prognostic value are only some of the aspects which caused the American Thoracic Society (ATS) to recommend this type of exercise testing.\*

### Smart solutions for an ideal test procedure

With Vyntus WALK the 6-minute walk test becomes truly mobile. To perform the test, only a pulse oximeter and a tablet PC are required. Real-time heart rate and  $\text{SpO}_2$  data are transferred and recorded via Bluetooth™. The manual entry of: blood pressure; oxygen supply; type of oxygen supply; as well as the rate of perceived exertion, allows for supplemental measurements and a plausibility check of your data. The tablet application has an intuitive workflow which smoothly moves you through the test procedure enabling you to focus on your patient. To ensure an optimal test process, handy pop-ups help you to guide and motivate the patient. The user-friendly interface lets you keep track of the measured values and the patient's stability. By means of the event button, events can be easily added during the test phase. As soon as the test is completed, they can be saved and transferred to the SentrySuite database together with the test data.

\* <https://www.thoracic.org/statements/resources/pfet/sixminute.pdf>



Promotional Material

**CPET Seminars and Workshops (Germany)**

Vyairé regularly offers different kinds of workshops and seminars:

**CPET**

Basic principles and measuring methods of resting and stress ECGs.

Cardiopulmonary Exercise Testing, device preparation and practical measurements.

Evaluation and simple interpretation.

**Practical Training**

Reliable settings, preparing and operating the measuring system.

Calibration, screen display, creating reports, load profiles, text editor, offline entries, predicted values as well as basic evaluation.

Cleaning the system and hygiene.

**Interpretation**

Principles of CPET. Assessment, evaluation and interpretation of measurements.

Ventilatory thresholds, pulmonary and/or cardiac limitations, blood gases and cardiac output as well as metabolism. Case studies.

Please ask for available appointments:

**CareFusion Germany 234 GmbH**

**Training Center**

Monika Schleyer  
Leibnizstrasse 7  
97204 Hoechberg  
+49 (0)931 4972-664 tel  
+49 (0)931 4972-382 fax  
training.eu@vyaire.com

Poster: Cardiopulmonary Exercise Testing; Typical normal CPET curves and curves in different diseases

Please call your local rep for assistance.

**CPET Training and Education (USA)**

Vyairé offers training and education at our Yorba Linda, CA training facility or on-site at your facility.

- Deliver comprehensive operational knowledge of the Vyntus™ CPX medical device and software.
- Emphasis placed on device understanding, measurement applications, American Thoracic Society (ATS) and American College of Sports Medicine (ACSM) clinical guidelines, customization, and quality control.
- Participants will be able to operate the Vyntus CPX device and system properly, safely, and effectively.

**Continuing Respiratory Care Education (CRCE) Credits**

- Many of our courses have been approved for CRCE credits by the American Association for Respiratory Care (AARC).

**Vyairé Medical**  
22745 Savi Ranch Parkway  
Yorba Linda, CA 92887  
714.919.3318 Tel

**Published by:**  
**CareFusion Germany 234 GmbH**  
Leibnizstrasse 7  
97204 Hoechberg  
+49 (0)931 4972-0 tel  
+49 (0)931 4972-423 fax  
www.vyairé.com

**Responsible:**  
Dr. Hermann Eschenbacher

Text and graphics have been compiled very carefully. Vyairé does not assume liability for mistakes and any damages resulting thereof. Subject to technical alterations.

Nothing in this publication is intended to be or substitute medical advice which may be provided only by qualified medical personnel.




The Canopy-option and the High / Low FIO<sub>2</sub>-option has not been and may not be cleared for commercial sale by the FDA Authority. Where applicable – country availability is dependent on the successful product registration with the National Authority of that country. Please read the complete Instructions For Use that come with the product.

All content and referenced material of this presentation is for informational purposes for the training participants only and are not intended to serve as a substitute for the consultation, diagnosis, and/or medical treatment by a qualified physician or healthcare provider.

## GLOBAL HEADQUARTERS

Vyairē Medical, Inc.  
26125 North Riverwoods Blvd  
Mettawa, IL 60045  
USA

 Vyairē Medical GmbH  
Leibnizstrasse 7  
97204 Hoechberg  
Germany



[vyairē.com](http://vyairē.com)

**Rx ONLY**

### For global distribution.

© 2019 Vyairē Medical, Inc. or one of its affiliates. All rights reserved. Vyairē and the Vyairē Medical logo are trademarks or registered trademarks of Vyairē Medical, Inc. or one of its affiliates. Medical devices class IIa (or correct class) according to Medical Devices Directive 93/42/EEC. Please read the complete Instructions For Use that come with the devices or follow the instructions on the product labelling.  
VYR-GLB-1900019