



Datt Mediproducts Private Limited



VELCARE[®] LYMPHEDEMA KIT

Bandage Set for Edema Management



PROPERTIES

- A compilation of all the necessary materials (including short stretch compression bandages, padding & fixation material)

MECHANISM

- Reduces the filtration rate
- Improves the efficiency of muscle & joint pumps
- Prevents the reaccumulation of evacuated lymph fluid
- Break up deposits of accumulated scar & connective tissue

CONTENT

Velnet™ (100% Cotton Stockinette)

- Used to protect the skin
- Made up of 100% Cotton which absorbs excess perspiration, it is important to use a stockinette underneath compression bandages

Digit Bandage (Woven Conforming Retention Bandage)

- Used to wrap individual digits
- Also used to keep foam rubber pads in place before short stretch bandages are applied

Velroll® Plus(orthopedic wool)

- Applied prior to foam Band
- Absorbs exudates & perspiration
- Evenly distribute pressure



INDICATION

- Treatment & Management of Lymphedema & associated disorders

FoamBand (Soft Foam Bandage)

- Applied just prior to the actual compression bandages
- To cushion the limb & prevent sharp indentations or irritation of the skin
- Serves to distribute the pressure evenly over area of arm or leg

• Velkomp (Short Stretch Compression Bandage)


- Used to achieve the required compression
- It's short stretch because of it's weave and it doesn't contain any elastomeric fibers

• Velcare™ Foam (High Density Foam Rubber)

- The use of pre-made or custom cut foam pads is recommended for anatomical indentations surrounding bony prominences (eg: Ankle & Wrist)
- Pieces of foam rubber are cut to size and placed in the indentation for even pressure distribution


Velpore™-SR (Surgical Paper Tape)

- Used to secure the short stretch bandages to avoid slippage for effective compression

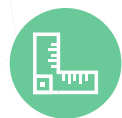

SIZES ARM KIT

Bandage	Sizes	Quantity
Digit Bandage	4 cm x 3.5 m	5
Velnet™	15 cm x 1.0 m	2
Velroll plus	10 cm x 3.0 m	1
FoamBand	10 cm x 2.5 m	1
VelKomp	6 cm x 5 m	1
VelKomp	8 cm x 5 m	2
VelKomp	10 cm x 5 m	2
Velpore™ -SR	1" x 5 m	2
Velcare Foam	9 cm x 1 cm	1
Carry Bag		1

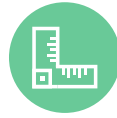

SIZES LEG KIT

Bandage	Sizes	Quantity
Digit Bandage	4 cm x 3.5 m	5
Velnet™	20 cm x 1.5 m	2
Velroll plus	15 cm x 3.0 m	1
FoamBand	10 cm x 2.5 m	2
VelKomp	8 cm x 5 m	1
VelKomp	10 cm x 5 m	3
VelKomp	12 cm x 5 m	2
Velpore™ -SR	1" x 5 m	2
Velcare Foam	9 cm x 1 cm	2
Carry Bag		1